



User Manual
6002009402
Getinge 86-series
S-8666, S-8668, S-8668T

The model name and serial number should always be provided when ordering spare parts or during queries per telephone or written correspondence:

Serial number of the machine:

.....

Machine model:

.....

The manual, or parts thereof, are a part of the product. This manual must not be modified, reproduced or copied without Getinge's written consent. The content may not be reused by third parties or used for unauthorized purposes. Violators will be prosecuted. Additional copies of this manual can be obtained from Getinge for a fee.

© Copyright 2013 Getinge. Reprinting prohibited.

Getinge Disinfection AB

Ljungadalsgatan 11

352 46 Växjö

SWEDEN

Internet: www.getinge.com

Contents

1	INTRODUCTION.....	7
1.1	About this manual.....	7
1.2	Terms used in this manual.....	7
1.3	Symbols in this manual.....	7
1.4	Intended use statement.....	7
1.4.1	Intended application areas.....	8
1.4.2	Intended user.....	8
1.4.3	Intended user environment.....	8
1.5	Potential misuse.....	9
1.6	Product liability.....	9
1.6.1	Declaration of Conformity.....	9
1.6.2	Disclaimer.....	9
1.6.3	Manufacturer.....	9
1.6.4	Warranty.....	9
1.7	Product labels.....	10
1.8	References.....	11
1.9	Service life.....	11
2	SAFETY.....	12
2.1	Safety functions.....	12
2.1.1	Door pinch protection.....	12
2.1.2	Password protection.....	12
2.1.3	Isolator switch.....	12
2.2	Safety statements.....	12
2.3	General warnings.....	13
2.4	Safety labels on the machine.....	13
2.5	Personal protection.....	13
2.6	Accident or incident reporting.....	14
3	DESCRIPTION.....	15
3.1	Configurations.....	15
3.2	Product description.....	16
3.2.1	Control panel.....	17
3.2.2	USB (option).....	19
3.2.3	Detergent cabinet.....	21
3.2.4	Guides for transport trolley.....	23
3.2.5	Lower guide for transport trolley.....	23
3.2.6	Outlet for external loading equipment (option).....	24
3.2.7	Guides for external loading equipment (option).....	24
3.2.8	Barcode scanner (option).....	24
3.2.9	Main power switch.....	24
3.2.10	Process stop (option).....	24
3.2.11	Door.....	24
3.2.12	Electrical outlet (option).....	25
3.2.13	Chamber.....	25
3.2.14	Carts.....	26
3.2.15	Automatic loading system (option).....	27

3.3	Control system.....	27
3.3.1	Access rights.....	28
3.3.2	Navigation and commands.....	29
3.3.3	Menu tree.....	31
3.3.4	Basic views.....	32
3.3.5	Change settings.....	35
3.3.6	View system information.....	37
3.3.7	Printing.....	37
3.3.8	System alerts.....	37
3.4	Accessories.....	40
3.5	Ambient illumination.....	40
3.6	Electromagnetic environment.....	40
4	OPERATION.....	41
4.1	Before using the machine.....	41
4.2	Workflow.....	41
4.3	Preparations.....	41
4.3.1	Loading the wash cart.....	42
4.4	Loading the chamber.....	43
4.4.1	Loading the chamber with an automatic loading system (option).....	43
4.5	Starting program.....	44
4.5.1	Starting program with a barcode scanner (option).....	44
4.5.2	Starting program with an automatic loading system (option).....	45
4.6	Unloading the chamber.....	45
4.6.1	Unloading the chamber with an automatic loading system (option)....	46
4.7	Check cleaning and disinfection results.....	46
4.7.1	Handling printouts (option).....	46
4.8	Holding or aborting program.....	47
4.9	Change program.....	48
4.10	Preparations for downtime.....	48
5	DETERGENTS.....	49
6	MAINTENANCE.....	51
6.1	Daily maintenance.....	51
6.2	Weekly maintenance.....	52
6.3	Monthly maintenance.....	52
6.4	Maintenance of parts inside the chamber.....	53
6.4.1	Maintenance of spray arms.....	53
6.4.2	Maintenance of spray arm manifold.....	54
6.4.3	Maintenance of strainer.....	55
6.5	Removing lime scale deposits in the chamber.....	56
6.6	Disinfecting chamber.....	57
6.7	Maintenance of the detergent system.....	57
6.7.1	Open the detergent cabinet door.....	57
6.7.2	Check for damage.....	57
6.7.3	Refill detergent.....	58
6.8	Maintenance of control panel and stainless steel surfaces	58
6.9	Maintenance of printer (option).....	59
7	TROUBLESHOOTING.....	60

7.1	System alerts.....	60
7.1.1	Handle alarms (login required).....	60
7.1.2	Handle warnings.....	60
7.1.3	Alarms and warnings list.....	61
7.1.4	Handle information message.....	68
7.2	General troubleshooting.....	70
7.3	Washing efficiency.....	71
7.4	Drying efficiency.....	71
7.5	Restart the machine.....	72
8	HANDLING OF WORN-OUT PRODUCTS.....	73
8.1	WEEE directive.....	73
8.2	Disposal of detergent or detergent containers.....	73
8.3	Disposal of the machine, accessories or spare parts.....	73
8.3.1	List of materials.....	74
9	ADDRESSES.....	75
	Alphabetical index.....	77

1 INTRODUCTION

This section contains information about the purpose and intended use of the machine. It also explains symbols and labels in the manual.

1.1 About this manual

This manual contains important information about usage of the machine. To ensure safe handling and error-free operation, read the manual before using the machine. Pay special attention to section 2 *SAFETY*. Some sections are not relevant for all machines; options, selections and customizations can vary.

1.2 Terms used in this manual

Name	Function/description
Load	Everything that is placed in the chamber, including wash carts and inserts.
Options	Custom choices that are not standard on the machine.
Soiled side	Loading side for soiled load.
Clean side	Unloading side for clean load.

1.3 Symbols in this manual

In this manual, warnings and information are indicated by the following symbols.

	Indicates a hazardous situation that could result in death, personal injury or damage to property if it is not avoided.
	Points out important details to consider for trouble-free and efficient usage of the machine.

1.4 Intended use statement

The intended use of the 86-series washer-disinfector is to clean, disinfect at a thermally low or intermediate level, and dry surgical instruments (solid and tubular), bowls, basins, glassware, receivers, suction bottles, baby bottles, anesthesia equipment and surgical shoes. The equipment can be made of stainless steel, aluminum, plastic, rubber, silicone or glass.

The device is not intended to clean and disinfect invasive devices as end point processing.

1.4.1 Intended application areas

The washer-disinfector performs the following basic functions:

- Cleaning of medical and laboratory equipment.
- Disinfection of medical and laboratory equipment.
- Drying of medical and laboratory equipment.

1.4.2 Intended user

The intended users of the 86-series washer-disinfector are operators and service personnel.

Operators

Operators are clinical personnel, using the machine for its intended use. Typical for this user group is a nurse or assisting nurse at a central sterile services department or personnel at an assistive device center or an animal facility for Life Sciences. The operator must understand the concepts of cleaning and disinfection as part of infection control in the intended use environment. The operator must have knowledge of local routines and requirements.

Typical tasks/responsibilities:

- Preparing, loading/unloading items, selecting and starting programs.
- Daily, weekly and monthly preventive maintenance to be performed by this category.

Service personnel

Typical for this user group is a Getinge technician or a local technician at the facility.

Typical tasks/responsibilities:

- Installation.
- Yearly preventive maintenance.
- Troubleshooting.
- Emergency shut downs.
- Education.
- End of life and recycling.

1.4.3 Intended user environment

The 86-series washer-disinfector is used in facilities such as: hospital wards, Central Sterile Supply Departments (CSSD), Theatre Sterile Supplier Units (TSSU) and hospital laboratories.

1.5 Potential misuse

- The washer-disinfector is not intended for use with household detergents.
- The washer-disinfector is not intended for cleaning and disinfection of bedpans or urine bottles (ISO 15883-3).
- The washer-disinfector is not intended for cleaning and chemical disinfection of thermolabile endoscopes (ISO 15883-4).

1.6 Product liability

1.6.1 Declaration of Conformity

The product compliance is stated in the Declaration of Conformity.

1.6.2 Disclaimer

- The manufacturer cannot guarantee that the product will work as intended if accessories or spare parts other than those approved by the manufacturer are used.
- The product liability only applies when the instructions in this manual have been followed in their entirety.
- The manufacturer accepts no responsibility for personal injuries when the instructions in this manual have not been followed in their entirety.
- The manufacturer reserves the right to change the specification and design of the product without any prior notice.

1.6.3 Manufacturer

This product was manufactured by:

Getinge Disinfection AB
Ljungadalsgatan 11
352 46 Växjö
SWEDEN




1.6.4 Warranty

The product warranty applies in accordance with the agreement reached between the Getinge sales company and the customer.

1.7 Product labels

The product labels are attached on the inside of the detergent cabinet. The product model and serial number are important when communicating with Getinge.

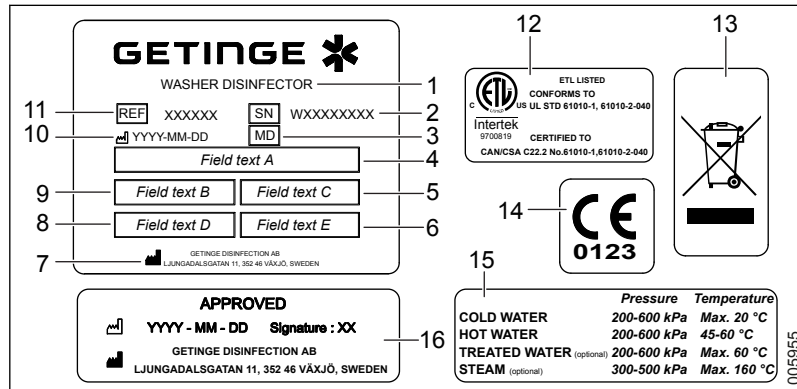

- | | |
|------------------------|-------------------------|
| 1. Designated purpose | 9. Maximum current |
| 2. Serial number | 10. Date of manufacture |
| 3. Medical device | 11. Machine model |
| 4. Voltage alternative | 12. Product listing |
| 5. Frequency | 13. Recycling symbol |
| 6. Protection class | 14. CE labeling |
| 7. Manufacturer | 15. Media connection |
| 8. Fuse | 16. Approval marking |

1.8 References

Documentation	Contents	Paper	Electronic media
User manual	Instructions for daily usage.	x	x
Installation manual	Instructions for assembly, installation and commissioning.	x	x
Quick guide	Instructions for frequently performed tasks.	x	x
Declaration of Conformity	Collection of ISO certificates.	x	-
Electrical diagrams	Collection of circuit diagrams.	-	x
Installation drawing	Installation measurements.	x	x
P&I diagrams	Collection of P&I diagrams.	x	x
Service manual	Instructions for service and preventive maintenance.	-	-
Spare parts	List of available spare parts.	-	-
Program sheets	Description of available programs.	-	x
Loading equipment for Getinge 86-series, 88 Turbo & CM320 washer-disinfectors.	Description of applicable loading equipment.	-	-
Smart loading trolley fixed height 2.0 for washer-disinfectors.	User manual for fixed height loading trolley (FHT Disinfection).	-	-
User manual Getinge loading equipment.	Instructions for daily use of Getinge loading equipment.	-	-
Loading guidebook for common hospital instruments in a Getinge washer-disinfector.	Description of how to load soiled surgical instruments for best possible cleaning results.	-	-

x = Delivered with the machine.

All documents are available at Getinge Extranet.

1.9 Service life

This washer-disinfector has an estimated service life of 10 years (40 000 cycles), provided that the recommended maintenance and service stated in the manuals have been carried out according to the instructions.

2 SAFETY

This section highlights important safety aspects that increases safety awareness. The user must follow the instructions and be familiar with the safety functions before using the machine.

2.1 Safety functions

2.1.1 Door pinch protection

If anything is detected between the door and the door frame when closing, the built-in door pinch protection activates and the door reopens. The door will not close as long as it is obstructed.

2.1.2 Password protection

The password protected programs can only be started by users with extended rights (login required). For more information, see section 3.3.1 *Access rights*.

2.1.3 Isolator switch

The machine must have a lockable isolator switch for the electric power supply (connected before/outside the machine's electrical connection). The isolator switch is installed by the customer and must be placed on a wall close to the machine. In an emergency, the main power can be switched off using the isolator switch.

2.2 Safety statements

Hazardous operations and situations are highlighted in the text.

Safety statement structure:	
• Warning symbol and signal word. • What is the risk? (Short description) • What causes the risk? (More precise information) • How can the risk be avoided?	
Example:	
	WARNING! Risk of burns <i>The inside surfaces of the chamber and the load may be hot directly after disinfection.</i> - Let the chamber cool down before entering the chamber. - Be careful when handling hot loads. - Wear protective clothing, such as gloves.

2.3 General warnings



WARNING!

Risk of electric shock and burns

There may be remaining voltage behind the cover panels, even when the power is switched off.

- Be careful when removing the cover panels.
- Installation and service behind the cover panels must only be performed by personnel with documented experience.



WARNING!

Risk of burns

The inside of the chamber and the load may be hot.

- Let the chamber and the load cool down before unloading.
- Wear protective gloves.



WARNING!

Risk of burns

Hot steam may escape from the chamber when the door is opened after a process is stopped or completed.

- Pay attention when standing next to the machine.

2.4 Safety labels on the machine

Symbol	Hazard	Measure
	Warning or Caution	Pay attention.
	Warning about hot surfaces	Pay attention when in contact with the washer-disinfector's surfaces.
	Warning about dangerous voltage	Pay attention when close to live parts. Do not touch electrical equipment or installations.

2.5 Personal protection

Working with the machine includes handling heavy loads and contact with hot surfaces and chemicals. Use protective gear, for example, footwear, gloves, clothing, glasses and other equipment according to local regulations.

2.6 Accident or incident reporting

Accidents or incidents must be reported according to local regulations.

3 DESCRIPTION


This section gives an overview of the machine and describes machine details and the control system. It contains necessary information for usage and understanding of the operating instructions.

3.1 Configurations

Model	Chamber width (mm)	Chamber length (mm)	Chamber height (mm)	Effective chamber volume (L)	Gross chamber volume (L)
S-8666	650	720	667	312	430
S-8668	650	810	667	351	480
S-8668T	650	810	667	351	480


3.2 Product description

The following illustration shows the machine as seen from the soiled side. Depending on the actual machine configuration, certain parts can appear on either the soiled side, the clean side or both.

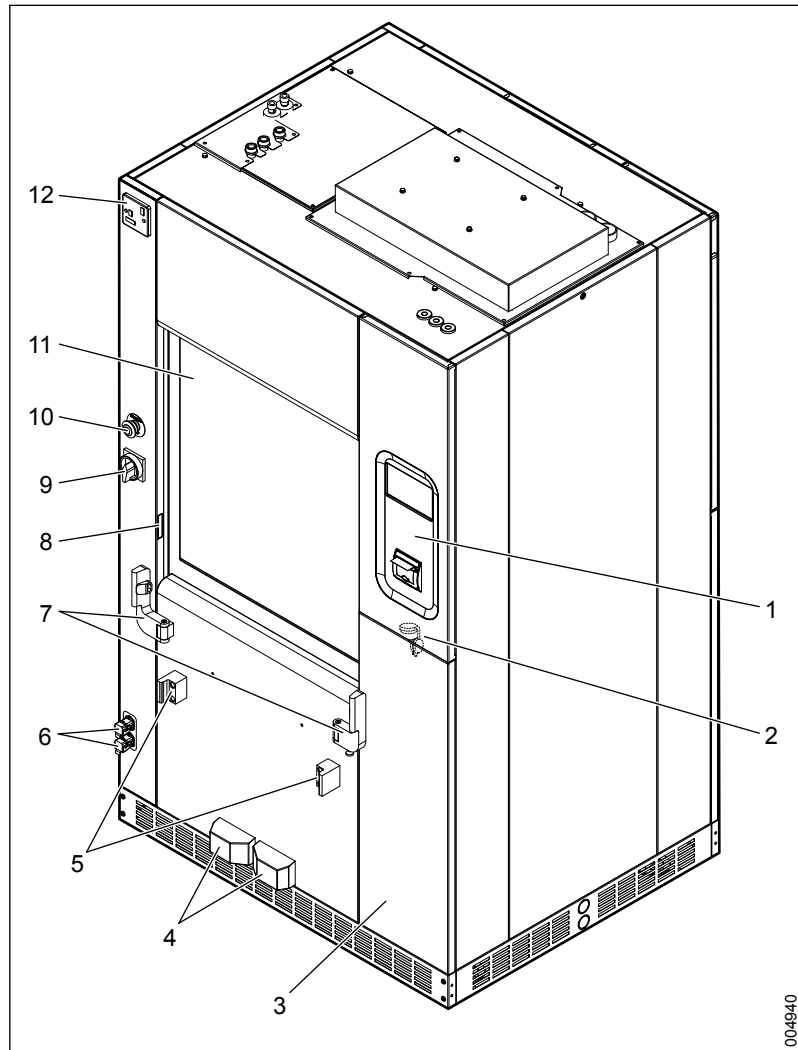

- | | |
|---|---------------------------------|
| 1. Control panel | 7. Guides for transport trolley |
| 2. USB (option), behind detergent cabinet door | 8. Barcode scanner (option) |
| 3. Detergent cabinet | 9. Main power switch |
| 4. Lower guide for transport trolley | 10. Process stop (option) |
| 5. Guides for external loading equipment (option) | 11. Door, chamber |
| 6. Outlet for external loading equipment (option) | 12. Electrical outlet (option) |

3.2.1 Control panel

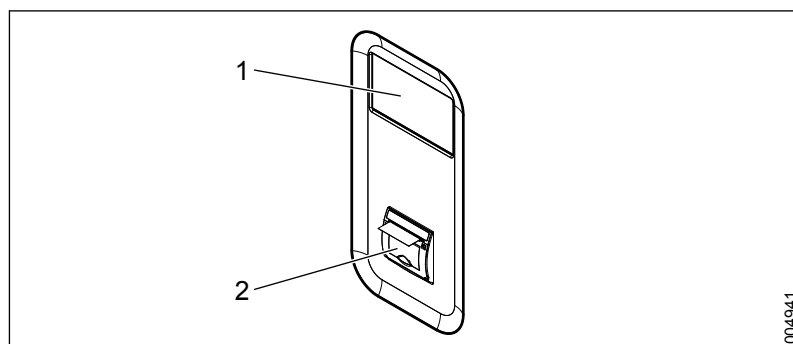

1. Touchscreen
2. Printer (option)




A machine with two doors has one control panel on each side (soiled side and clean side). The control panel layout may vary depending on the options installed on the washer-disinfector.

The machine is controlled and monitored from a control panel with a pressure-sensitive touchscreen. For more information, see section 3.3 *Control system*.

3.2.1.1 Touchscreen

The machine is controlled by tapping commands on the touchscreen. The touchscreen also displays alerts. There is one touchscreen on each control panel. For more information, see section 3.3 *Control system*.

3.2.1.2 Printer (option)

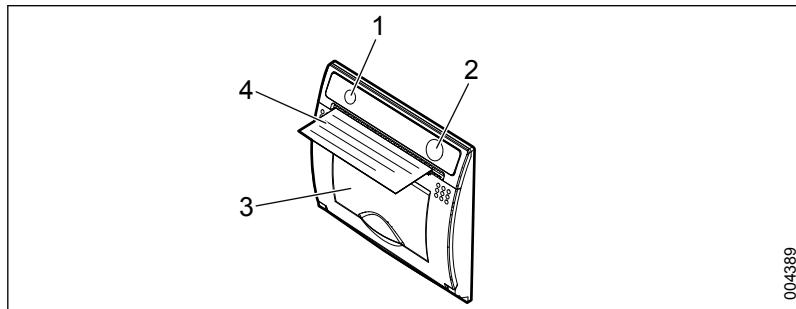

- | | |
|----------------------|----------|
| 1. Indication lamp | 3. Lid |
| 2. Feed paper button | 4. Paper |


After a program, the printer will automatically print reports (log data). In case of machines with double doors, the printer is only installed on one control panel, either on the clean or on the soiled side.

The indication lamp shows the status of the printer:


- A steady green light = the printer is working correctly.
- A flashing green light = there is something wrong with the printer, for example, the paper needs to be refilled or aligned.

For more information, see section 3.3.7 *Printing* and 6.9 *Maintenance of printer (option)*.

3.2.2 USB (option)

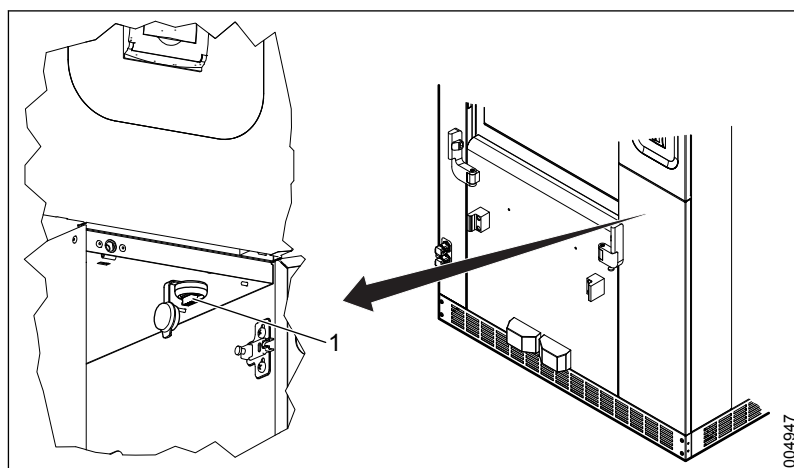

1. USB

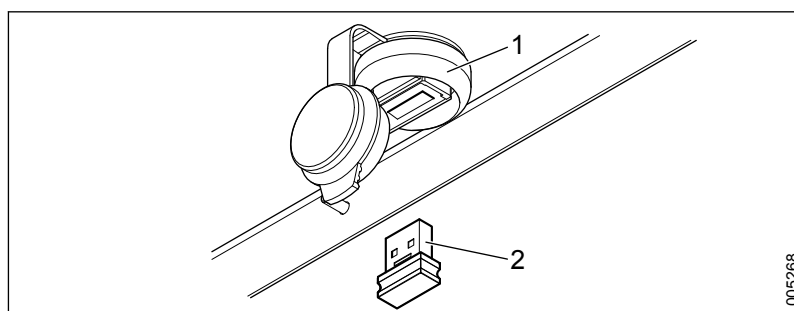

1. USB socket
2. USB stick



Only use the USB stick delivered by Getinge.




Do not remove the USB stick while data is being saved. When the touchscreen displays a warning that the USB memory is almost full, it is time to transfer the data to another storage device.

The USB option is used as a backup for process reports. The USB port is located inside the detergent cabinet. For more information, see section 6.7.1 *Open the detergent cabinet door*.

The activation of the USB function and the USB stick is optional. If activated, the process reports are saved automatically on the USB stick. A flashing light on the USB stick indicates that data is being saved.

When the touchscreen displays a warning that the USB memory is almost full, it is time to transfer the data to another storage device.

3.2.3 Detergent cabinet

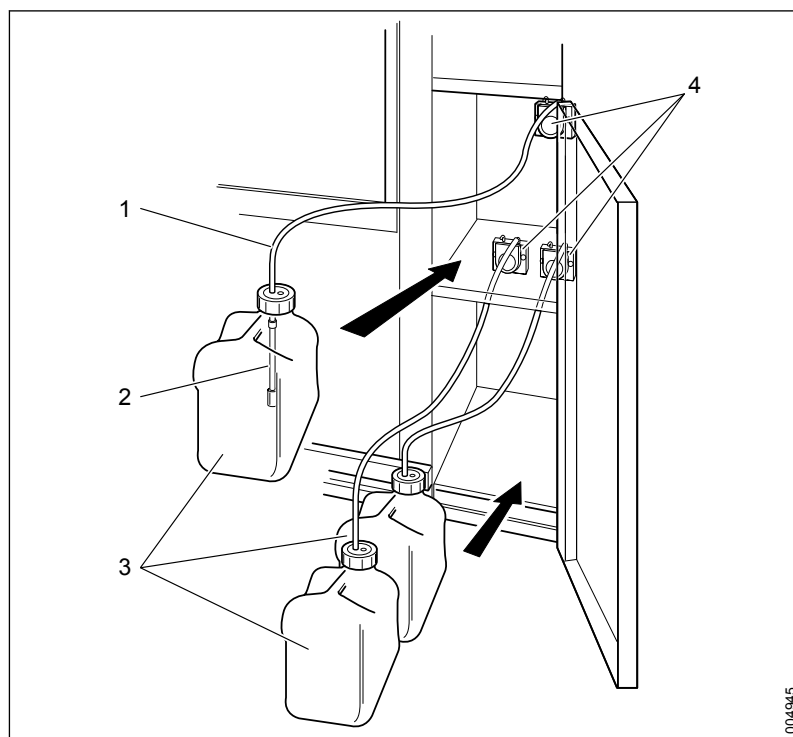

- | | |
|-------------------|-------------------------|
| 1. Transport hose | 3. Detergent containers |
| 2. Suction lance | 4. Dosing pumps |

The detergent cabinet contains the detergent system that doses and distributes the detergent to the chamber. As standard, two dosing pumps are installed. Two extra dosing pumps are optional.

As standard, the system issues a warning if the detergent level is too low in any of the containers.

As an option, flow monitoring can be installed. This system issues a warning if the detergent flow between the container and the chamber is too low.

For more information, see section 6.7.1 *Open the detergent cabinet door*.


3.2.3.1 How detergents are labeled



If a detergent label or dosage tag is missing or if the product name on the label is not the same as the product name on the detergent container, contact service personnel.

Using Getinge detergents


The detergent label is attached to the suction lance. The label corresponds to the name on the detergent container. The dosage tag is attached to the cable close to the detergent label.

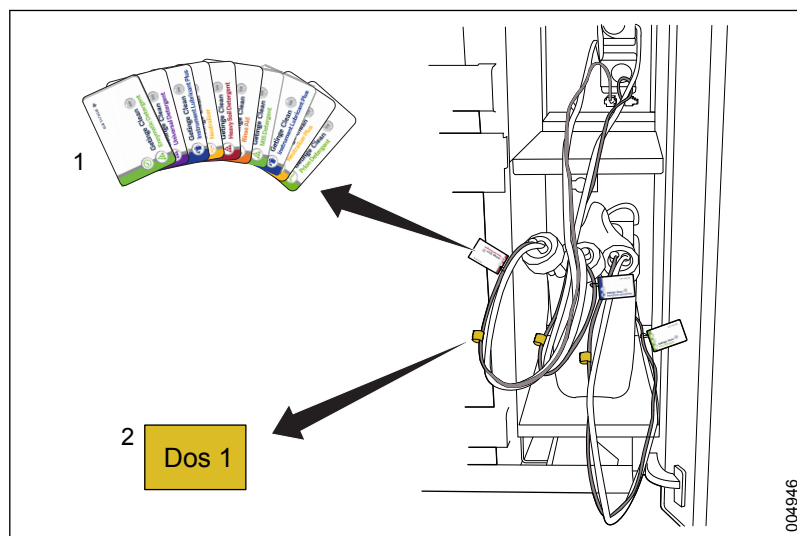

1. Detergent label
2. Dosage tag

Using non-Getinge detergents


Non-Getinge detergents do not have a detergent label. The detergents used are listed with names on the list of detergents. The names on the list correspond to the names on the detergent containers. The dosage tag is attached to the cable close to the detergent label.

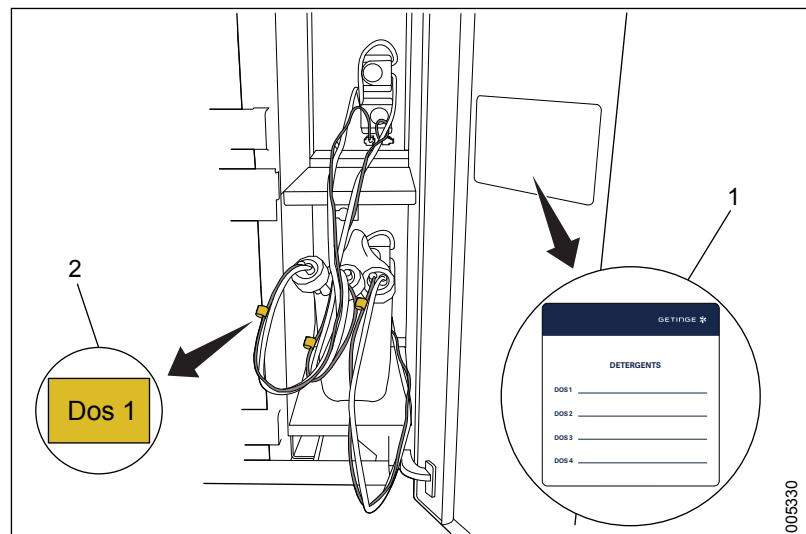

1. List of detergents
2. Dosage tag

3.2.4 Guides for transport trolley

The guides ensure that the transport trolley is in the proper position on the loading side. Guides cannot be combined with an external loading system or AGS.

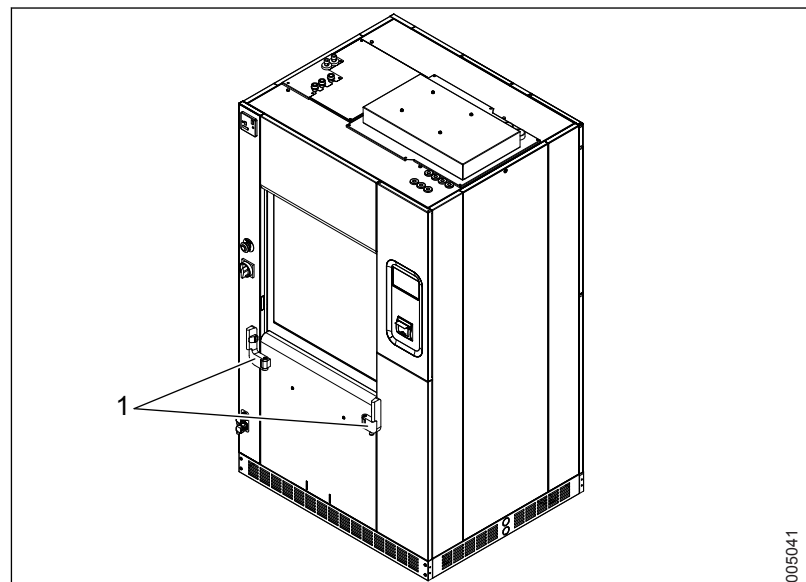

1. Guides for transport trolley

3.2.5 Lower guide for transport trolley

The lower guide for transport trolley can only be used in combination with external loading equipment. It is used to ensure correct positioning when loading manually.


3.2.6 Outlet for external loading equipment (option)

The outlet is used to connect external loading equipment.

3.2.7 Guides for external loading equipment (option)

The guides ensure correct positioning of the external loading equipment when loading and unloading.

3.2.8 Barcode scanner (option)

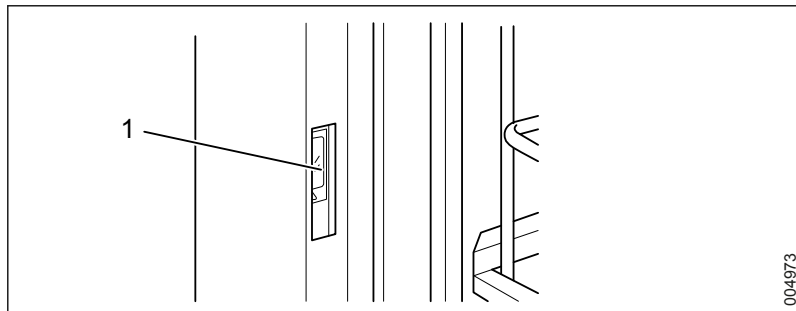

1. Barcode scanner

The built-in barcode scanner (option) is used for automatic program selection. The barcode scanner reads the barcode on the wash cart as it enters the chamber.

3.2.9 Main power switch

The main power switch is used to turn on or off the device. For more information, see section 7.5 *Restart the machine*.

3.2.10 Process stop (option)

The process stop button is an option only available for the UK market. It is used to stop the program (without switching the power off).

Activating the process stop triggers an alarm that must be dealt with as indicated in section 3.3.8.3 *Alarms*.

3.2.11 Door

The machine can be equipped with one or two doors:

- One door for loading and unloading on the same side.
- Two doors, one door for loading and one door for unloading on the opposite side.

The doors are motorized and are opened and closed vertically.

The doors have a door pinch protection. For more information, see section 2.1.1 *Door pinch protection*.

On machines with two doors, the interlocking system prevents cross-contamination from the soiled side to the clean side by only allowing one door to open at a time.


3.2.12 Electrical outlet (option)

The electrical outlet is an option only available for the UK market. Residual current device included.

Voltage (V)	Frequency (Hz)	Current (A)	Residual current device (mA)
230	50	10	30

3.2.13 Chamber

The chamber volume differs between different machine models. For more information, see section 3.1 *Configurations*.

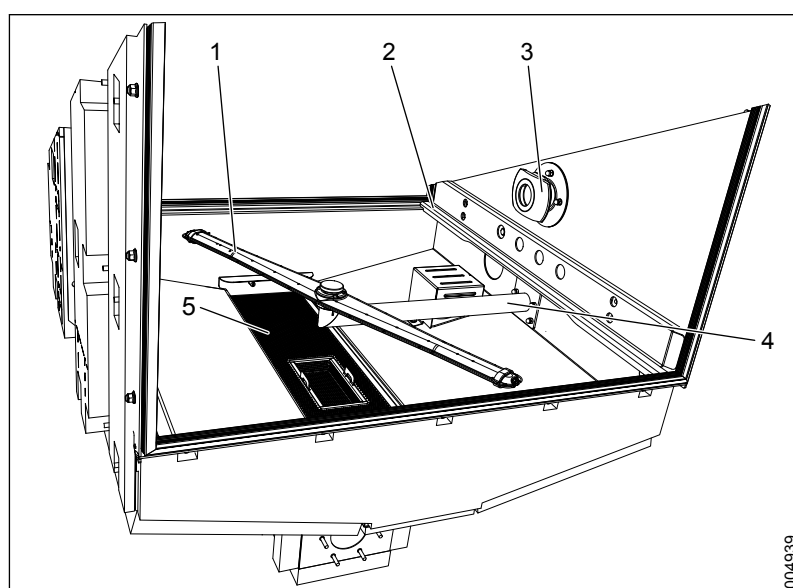

- | | |
|---------------|-----------------------|
| 1. Spray arm | 4. Spray arm manifold |
| 2. Guide rail | 5. Strainer |
| 3. Docking | |

3.2.13.1 Spray arm

During the cleaning phase, the spray arm rotates and distributes water and detergent evenly in the chamber. During the drying phase, the spray arm distributes heated air.

When SWRSS (option) is installed, the spray arm rotation is monitored, and the system warns for deviations.

3.2.13.2 Guide rail

The horizontal metal tracks guide the wash cart into the chamber during loading.

3.2.13.3 Docking

The docking secures and keeps the wash cart in place in the chamber and acts as an inlet for water and heated air.

3 DESCRIPTION

3.2.13.4 Spray arm manifold


The spray arm manifold distributes water and detergent or heated air to the spray arm.

3.2.13.5 Strainer

The strainer filters large particles in the water to allow the water to circulate and to avoid conditions for growth of bacteria.

3.2.14 Carts

Wash carts and transport trolleys are not included in the machine configuration, but are necessary for a correct usage of the machine.

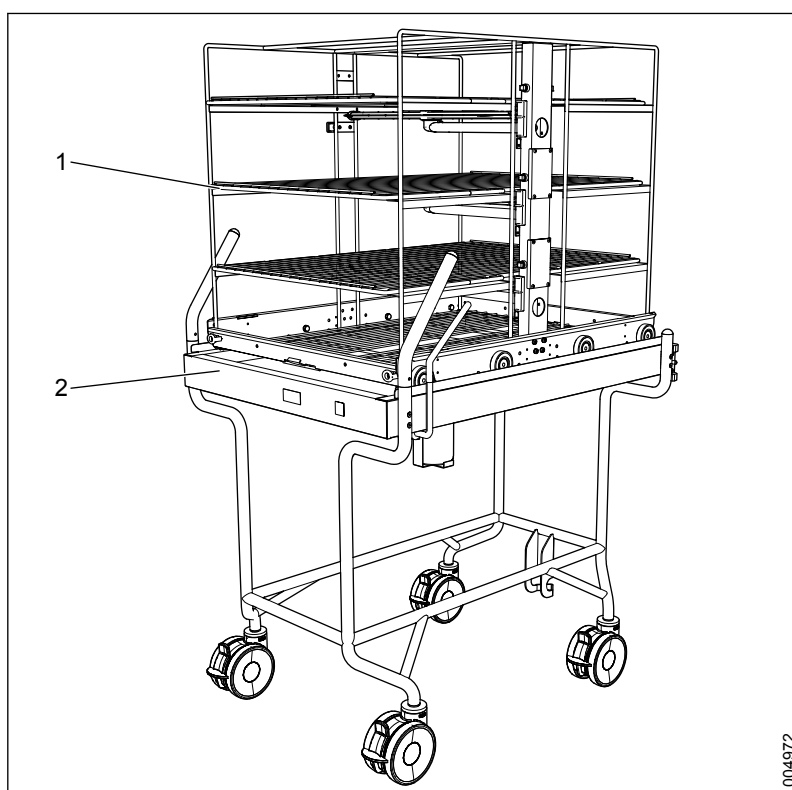

1. Wash cart
2. Transport trolley

3.2.14.1 Wash cart

Place the soiled items in the wash cart and load the wash cart into the chamber.

There are various types of wash carts, depending on the type of soiled load. For more information about the different loading systems, see the manual *Loading equipment for Getinge 86-series, 88 Turbo & CM320 washer-disinfectors*.

3.2.14.2 Transport trolley

The transport trolley is used to transport the wash cart.

The transport trolley has a docking function that locks the transport trolley to the docked unit when in the correct position.

There is a handle on the transport trolley to release it from the docked unit.

For more information, see the manual *Smart loading trolley fixed height 2.0 for washer-disinfectors*.

3.2.15 Automatic loading system (option)

The washer-disinfector can be configured for manual or automated loading/unloading, or a combination of both.

Machines intended for the automatic loading system are fitted with a fixed scanner for automatic program selection.

Free-standing loading and unloading conveyors (FSLC, FSUC) are part of an automatic system for loading/unloading of wash carts. Getinge AGS 2.0 is a fully automatic and highly flexible loading and unloading system. For more information about the loading systems, contact your local Getinge sales representative.

3.3 Control system

The user interface consists of a control panel with a touchscreen. Programs are started and settings are adjusted by tapping touch buttons on the touchscreen. The machine can be equipped with one or two control panels. For more information, see section 3.2.1 *Control panel*.

3.3.1 Access rights

The daily operator access rights are referred to as default and do not require a login.

Default access rights are granted the following basic functions:

- Selecting programs
- Opening/closing doors
- Starting programs
- Confirming warnings.

The extended access rights are referred to as user and require a login.

User rights grant access to the following additional functions:

- Running password-protected programs
- Confirming alarms
- Printing emergency printouts.



Factory settings:

User name = User

Password = 1234


The factory settings may be changed during commissioning.

3.3.2 Navigation and commands

3.3.2.1 Keyboard

The display allows the operator to type letters, special characters and numeric values.

The virtual keyboard appears automatically when the text field is tapped.

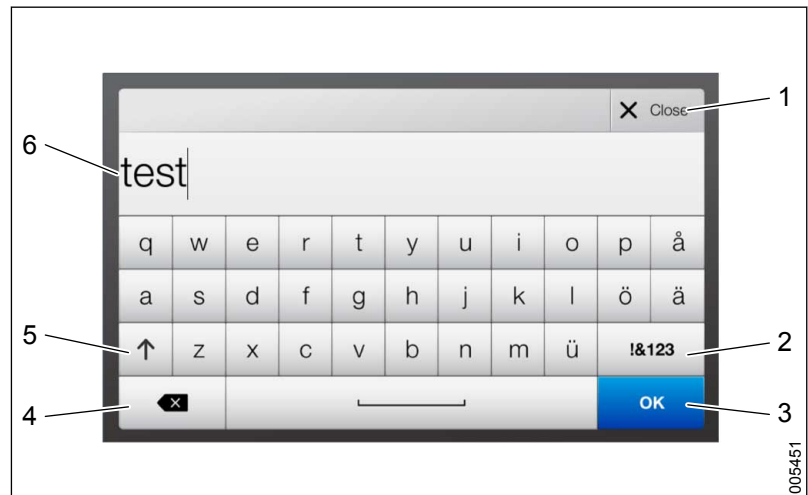

- | | |
|-----------------------------------|------------------------|
| 1. Close | 4. Erase text |
| 2. Numeric or alphabetic keyboard | 5. Upper or lower case |
| 3. OK | 6. Text field |

Name	Description
Close	Button for closing the window and returning to the previous view.
Numeric or alphabetic keyboard	Button for switching between numeric and alphabetic keyboard.
OK	Button for confirming the new setting.
Erase text	Button for erasing text in the text field.
Upper or lower case	Button for switching between upper and lower case letters in the text field.
Text field	Field for inserting text.


3.3.2.2 Touchscreen commands

Command	Symbol	Usage/function
Abort		Abort process
Back		Return to previous window
Close		Close current window.
Close Door		Close the door
Erase		Erase text inserted with keyboard
Information		View detailed information
Menu levels	● = level one ○● = level two ○○● = level three	Show current level of submenu
OK		Confirm change
Open Door		Open the door
Save		Save changes
Scroll down		Scroll down in lists
Scroll up		Scroll up in lists
Start Program		Start a selected program

3.3.3 Menu tree



The available menus depend on the user's access rights.

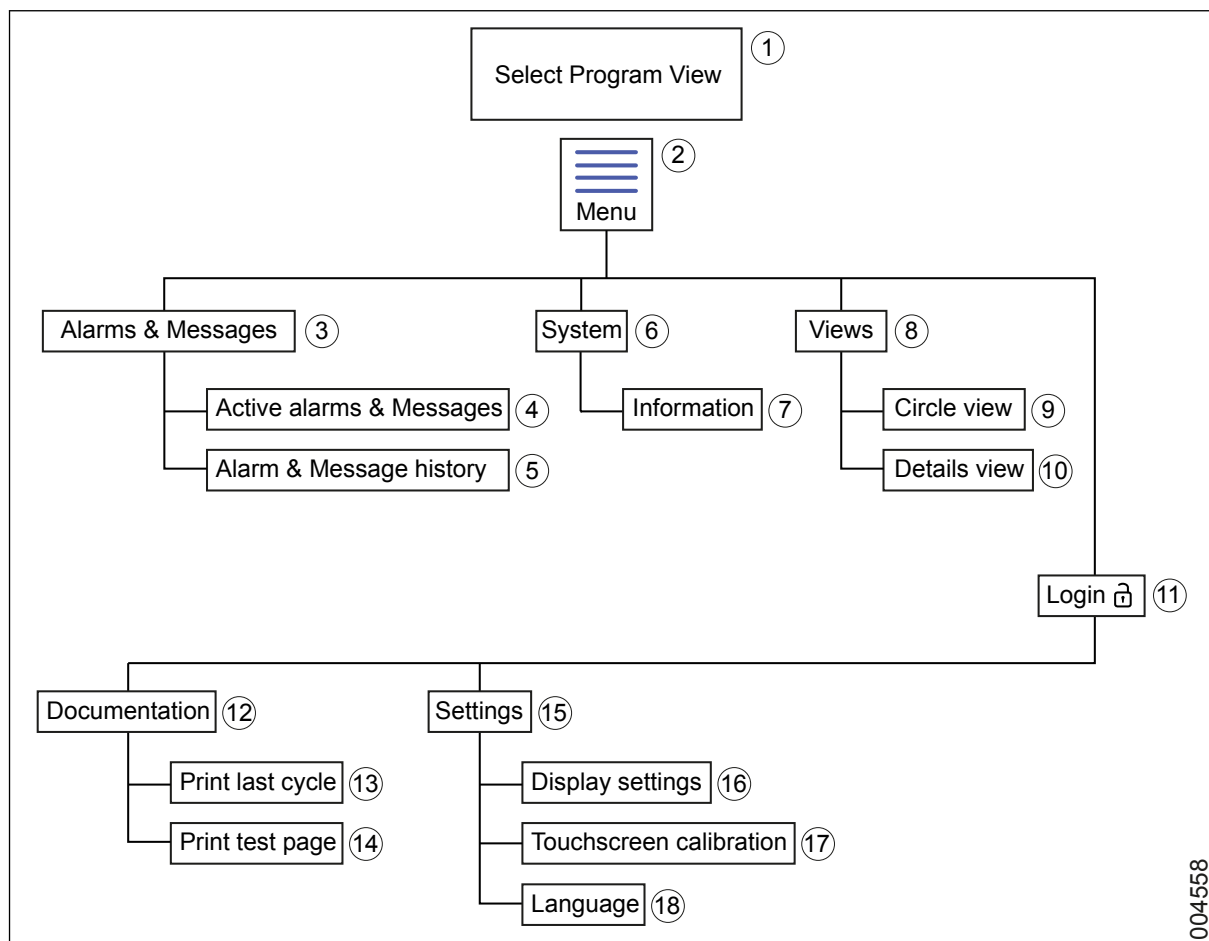

- | | |
|---|------------------------------------|
| 1. Select Program view | 10. Details View |
| 2. Menu | 11. Login |
| 3. Alarms & Messages | 12. Documentation |
| 4. Active Alarms & Messages | 13. Print Last Cycle |
| 5. Alarms & Messages History | 14. Print Test Page |
| 6. System | 15. Settings |
| 7. Information | 16. Display Settings |
| 8. Views | 17. Touchscreen Calibration |
| 9. Circle View | 18. Language |


004558

3.3.4 Basic views

The buttons on the touchscreen are always visible, but are only available for selection when highlighted. For example, the **Open Door** button is not highlighted during a program since it is not possible to open the doors.

3.3.4.1 Select Program view (start view)

When the machine is ready for use, the **Select Program** view (start view) is displayed on the touchscreen.

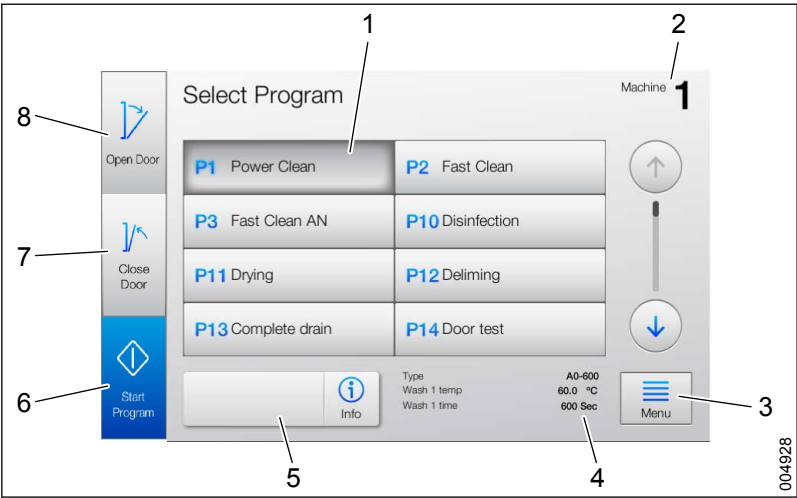

- | | |
|-----------------|------------------|
| 1. Program list | 5. Info |
| 2. Machine | 6. Start program |
| 3. Menu | 7. Close Door |
| 4. Quick view | 8. Open Door |


Name	Description
Program list	Buttons for selecting a program.
Machine	Displays the machine ID.
Menu	Button for displaying sub menus.
Quick view	Displays a selection of program parameters.
Info	Button for displaying detailed information.
Start program	Button for starting a program.
Close Door	Button for closing the door.
Open Door	Button for opening the door.

3.3.4.2 Menu

Path: Select Program > Menu



The available menus depend on the user's access rights.


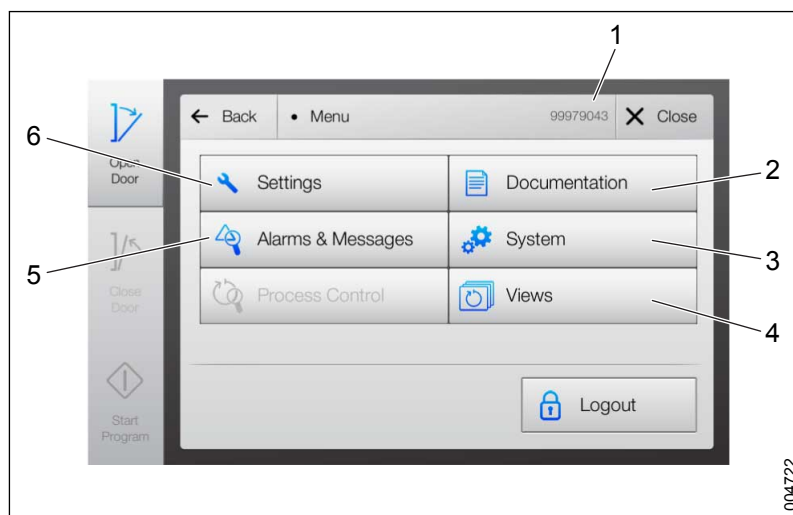


- | | |
|------------------|----------------------|
| 1. User | 4. Views |
| 2. Documentation | 5. Alarms & Messages |
| 3. System | 6. Settings |

Name	Description
User	Displays the currently logged on user.
Documentation	Button for printing reports.
System	Button for displaying system information.
Views	Button for Status View settings.
Alarms & Messages	Button for displaying alerts and warnings information.
Settings	Button for touchscreen settings.

3.3.4.3 Status view (during program)

The **Status View** is displayed during a currently running program. It contains information about the progress. The **Status View** can be displayed in **Circle View** or **Details View**.

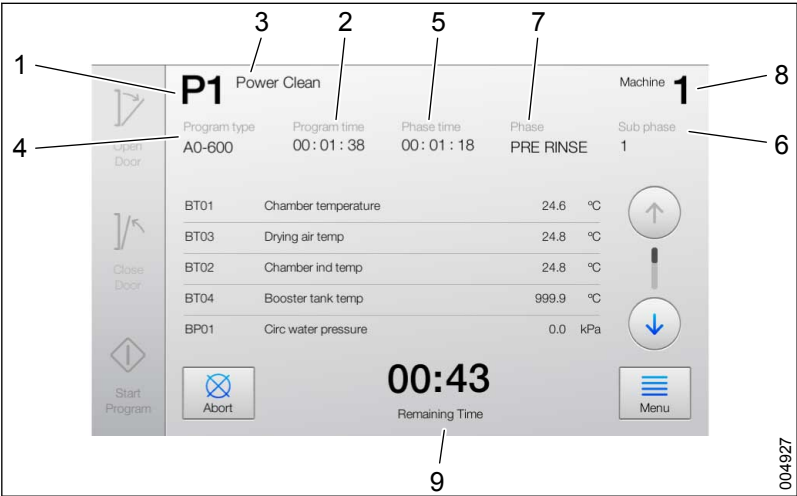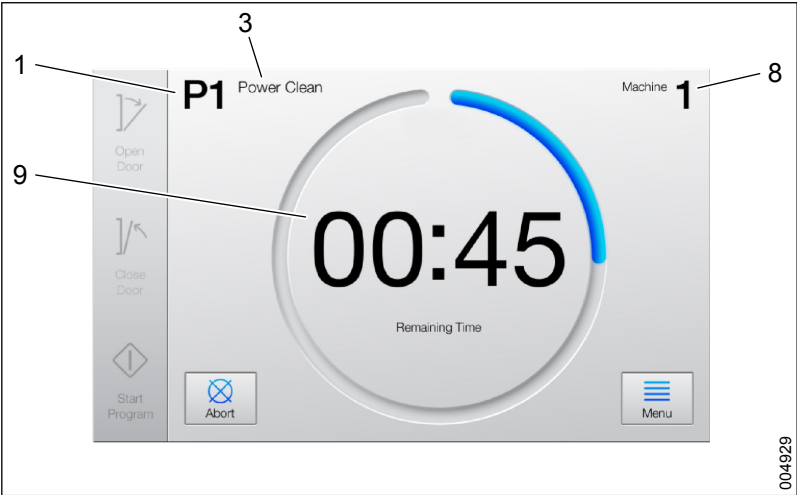

- | | |
|-------------------|-------------------|
| 1. Program number | 6. Sub phase |
| 2. Program time | 7. Phase |
| 3. Program name | 8. Machine |
| 4. Program type | 9. Remaining Time |
| 5. Phase time | |

Name	Description
Program number	Displays the program number.
Program time	Displays the total program time.
Program name	Displays the program name.
Program type	Displays the program type.

Name	Description
Phase time	Displays the elapsed time of the current phase.
Sub phase	Displays the current sub-phase number.
Phase	Displays the phase name.
Machine	Displays the machine ID.
Remaining Time	Displays the remaining program time.

3.3.5 Change settings

Some settings require login. For more information, see section 3.3.1 *Access rights*.

3.3.5.1 Set status view

The default setting is **Circle View**.

Switch to a detailed view:

Path: **Menu > Views > Details View**

1. In the **Status View** tap **Menu** to change to **Details View**.

To go back without changing, tap **Abort**.

3.3.5.2 Set screensaver (login required)

Path: **Menu > Settings > Display Settings**

1. Change the time for the screensaver by tapping **minutes** field and enter a new number.
2. Tap **Save** to confirm and go back.


3.3.5.3 Set touchscreen brightness (login required)

Path: **Menu > Settings > Display Settings**

1. Adjust the display brightness by tapping **+** or **-**.
2. Tap **Save** to confirm and go back.

3.3.5.4 Calibrate touchscreen (login required)

A calibration of the touchscreen is necessary when the touchscreen is not responding correctly when tapping.

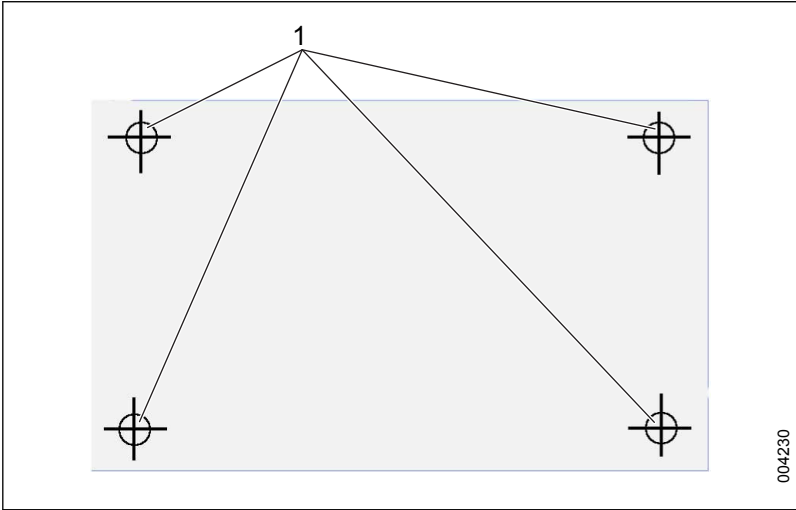

1. Crosshairs

Name	Description
Crosshairs	Tap to calibrate the touchscreen.

Path: **Menu > Settings > Touchscreen Calibration**

1. Tap the crosshairs, one at a time, as they appear.

When all four crosshairs have been tapped the calibration is saved automatically. If the calibration fails, it will not be saved and the previous settings are kept.

3.3.5.5 Set touchscreen and print-outs language (login required)

Path: **Menu > Settings > Language**

1. Tap the desired language in the list. Available languages are determined by the current machine configuration.
2. Tap **Save**.

3.3.6 View system information

3.3.6.1 View alarm and message information

The control system displays information about alarms, such as the name, alarm code, date, and time.

View currently active alarms:

Path: **Menu > Alarms & Messages > Active Alarms & Messages**

View old (not active) alarms:

Path: **Menu > Alarms & Messages > Alarms & Messages History**

3.3.6.2 View system and machine information

Displays important system parameters to be used by service personnel.

Path: **Menu > System > Information**

3.3.7 Printing

3.3.7.1 Print emergency print-outs (login required)

It is possible to print a report from the last program if, for example, the printer has run out of paper during a program.

Path: **Menu > Documentation > Print Last Cycle.**

1. Tap **Print**.

3.3.7.2 Print test page (login required)

Path: **Menu > Documentation > Print Test Page.**

1. Tap **Print**.


3.3.8 System alerts

The control system sends alerts with an information message, a warning or an alarm.

- An information message provides guidance before starting.
- A warning provides information about minor problems during a program or during standby.
- An alarm provides information about serious problems during a program or during standby. It also generates an alarm signal.

3.3.8.1 Information text

Information messages display conditions that must be resolved before starting a program.

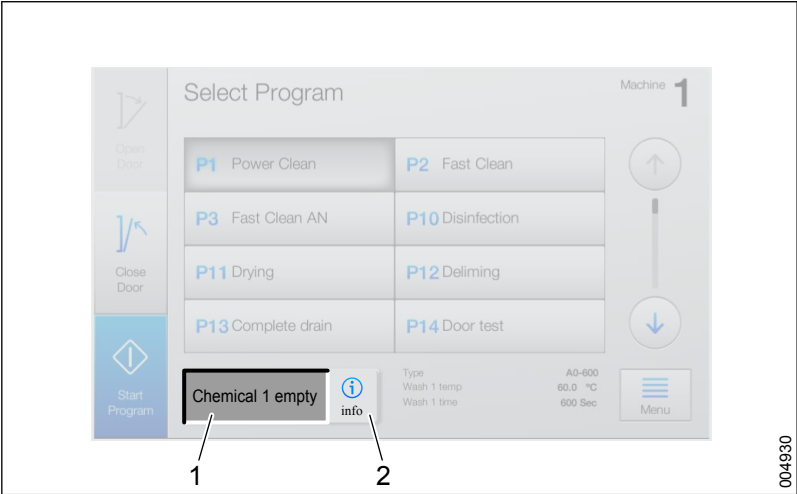


- 1. Information message
- 2. Info

Name	Description
Information message	Displays information about starting issues.
Info	Button for displaying detailed information.

For more information, see section 7 *TROUBLESHOOTING*.

3.3.8.2 Warnings

When a minor problem occurs, a yellow warning message appears on the touchscreen. The program continues.

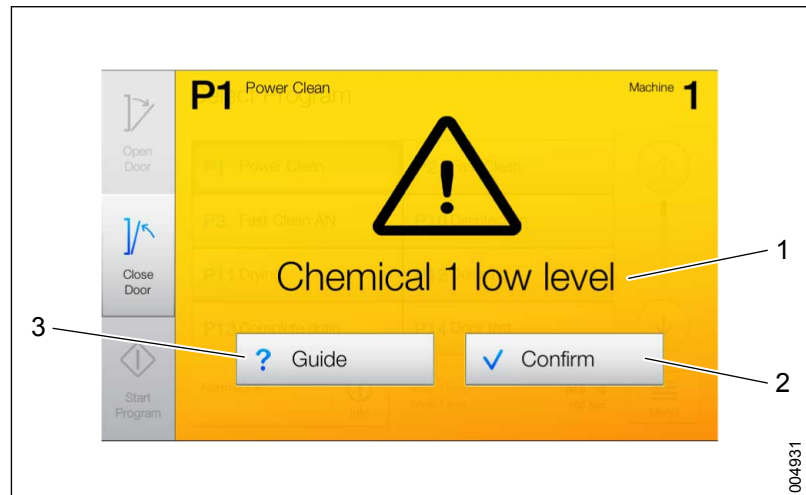


- 1. Warning message
- 2. Confirm
- 3. Guide

Name	Description
Warning message	Displays a warning message.
Confirm	Button for confirming the warning.
Guide	Button for displaying detailed information and guidance instructions.

For more information, see section 7 *TROUBLESHOOTING*.

3.3.8.3 Alarms

When a serious problem occurs, a red alarm message appears on the touchscreen. The program is stopped and an alarm signal sounds.

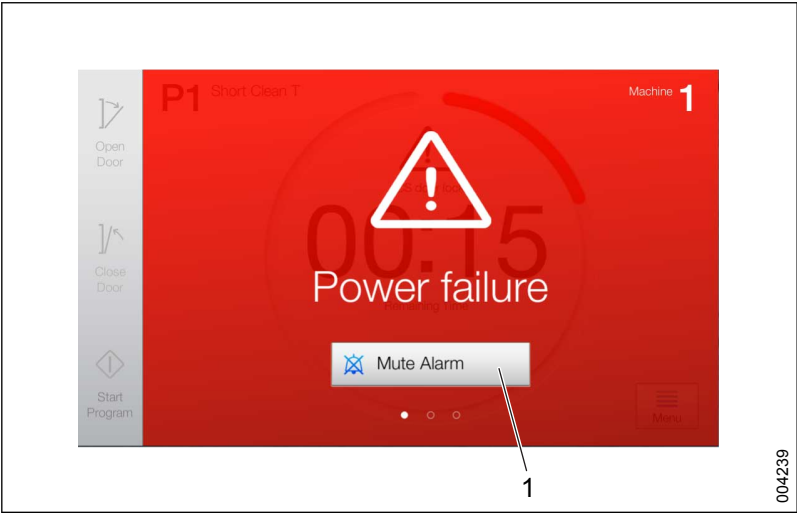

1. Mute Alarm

Name	Description
Mute Alarm	Button for muting the alarm signal.

For more information, see section 7 *TROUBLESHOOTING*.

3.4 Accessories

There are several accessories intended for use with this machine, such as wash carts. For a complete list, see the *Getinge loading equipment* manual.

3.5 Ambient illumination

The ambient illumination must be set according to ISO 12100.

3.6 Electromagnetic environment

The washer-disinfectors in the Getinge 86-series comply with the EMC standard.

4 OPERATION

This section contains operating instructions and routines for daily usage of the machine.

4.1 Before using the machine

If the machine has not been used for more than 72 hours, the program **Self-disinfect** should be run. It is recommended to run a program with an empty chamber before using the machine to clean items.

4.2 Workflow

The following steps describes the basic workflow:

1. **Preparations**
Prepare and pre-clean the items. Load the wash cart. For more information, see section 4.3 *Preparations*.
2. **Loading the chamber**
Load the chamber with or without an automatic loading system (option). For more information, see section 4.4 *Loading the chamber*.
3. **Starting program**
Select and start a program manually, with a barcode scanner (option), or with an automatic loading system (option). For more information, see section 4.5 *Starting program*.
4. **Unloading the chamber**
Unload the chamber with or without an automatic loading system (option). For more information, see section 4.6 *Unloading the chamber*.
5. **Check cleaning and disinfection results**
Inspect the items. For more information, see section 4.7 *Check cleaning and disinfection results*.

4.3 Preparations

1. Remove non-washable parts from the load.
2. Pre-clean the items. Use pre-treatment agents if the items are very soiled or have been left to dry for several hours.
3. Ensure that no detergent or pre-treatment agent is left from pre-cleaning.

4.3.1 Loading the wash cart



WARNING!

Risk of contamination due to cut injuries

Exposure to infected loads could lead to life-threatening conditions.

- Place sharp and pointed utensils in a way that they cannot cause injuries when loading and unloading.

Load the wash cart and observe the following guidelines:

- Always follow the cleaning instructions from the manufacturer of items and instruments.
- Empty all containers that are to be washed before they are put in the wash cart.
- Only use inserts intended for the items and utensils.
- Place sharp and pointed utensils in a way that they cannot cause harm when they are loaded and unloaded.
- Ensure that the utensils are not placed on top of each other, preventing the water from reaching all surfaces.
- Place containers with openings downwards in the intended inserts. The water must be able to flow freely through the vessels.
- Tilt items in a way that prevents collection of water.
- Rigid endoscopes must be dismantled in accordance with the manufacturer's instructions before cleaning.
- Open all valves and remove covers and seals on items in a way that the channels can be thoroughly flushed out.
- Check that the items are retained in their inserts.
- Distribute the load evenly.
- Ensure that the water can flow through perforated instruments. Special accessories may be needed for cleaning certain instruments.
- Ensure that the spray arms can rotate freely without touching any of the items.
- Ensure that the spray arms reach all surfaces of the items.
- A fine filter must be used to clean and disinfect tubular instruments with a passage of less than 2 mm. The fine filter is located at the inlet to the wash ramp.

For more information about how to load the wash cart, see the manual *Loading guidebook for common hospital instruments in a Getinge washer-disinfector*.

4.4 Loading the chamber

WARNING!


Risk of contaminated items



Incorrect loading of items or improper use of accessories can affect the cleaning result.

- Only use approved accessories.
- Follow manufacturer's reprocessing instructions for medical items.

1. Open the door on the soiled side by tapping **Open Door**.
2. Before loading, ensure that there are no items on the chamber bottom.
3. Dock the transport trolley. For more information, see section 3.2.14.2 *Transport trolley*.
4. Push the wash cart into the chamber.

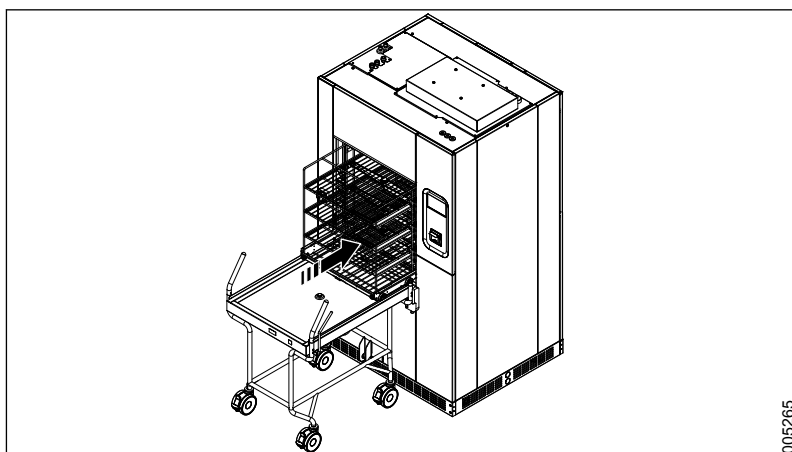

5. Undock the transport trolley.

4.4.1 Loading the chamber with an automatic loading system (option)

1. Dock the transport trolley. For more information, see section 3.2.14.2 *Transport trolley*.
2. Push the wash cart onto the loading station. The loading station starts and actively moves the wash cart to its holding position or into the chamber.
3. Undock the transport trolley.

4.5 Starting program



WARNING!

Risk of contaminated items

Incorrect selection of a program can affect the cleaning result.

- Follow operation instructions in User manual.
- Visually inspect the items after unloading.
- Follow manufacturer's reprocessing instructions for medical items.



During commissioning, changes may have been carried out to default program names, parameter settings, password protection usage and other settings.



Program names and order of appearance may be different from the pre-defined names, and the range of program may be limited or expanded.

The machine is equipped with factory installed standard programs. For more information, see the Program sheets.

1. Select a program on the touchscreen.
2. Tap **Start program** to close the door and to start the selected program automatically.

Certain programs require a login for access. For more information, see section 3.3.1 *Access rights*.

4.5.1 Starting program with a barcode scanner (option)

The machine may be equipped with a barcode scanner (option) for automatic program selection. As the wash cart enters the chamber, the barcode scanner reads the barcode on the wash cart and selects the correct program.

1. Tap **Start program**.
The door closes and the program status is displayed on the touchscreen.

It is possible to manually select and run a different program, other than the one that is automatically read via the barcode scanner.

1. Select the program to be used instead of the automatically selected.
2. Tap **Start program**.

4.5.2 Starting program with an automatic loading system (option)

The barcode scanner automatically selects a program by reading the barcode on the wash cart as the cart enters the chamber. When a wash cart is loaded into the chamber, the door closes and the selected program starts automatically. The program status is displayed on the touchscreen.

4.6 Unloading the chamber

WARNING!

Risk of burns




The surface on items and accessories may be hot.

- Let the items and accessories cool down before handling.
- Wear protective gloves.

When the program is completed, the touchscreen displays **Complete** and the door opens.

1. If the touchscreen has gone into idle mode, tap on it to reactivate.
2. Dock the transport trolley. For more information, see section 3.2.14.2 *Transport trolley*.
3. Unload the wash cart from the chamber onto the transport trolley.

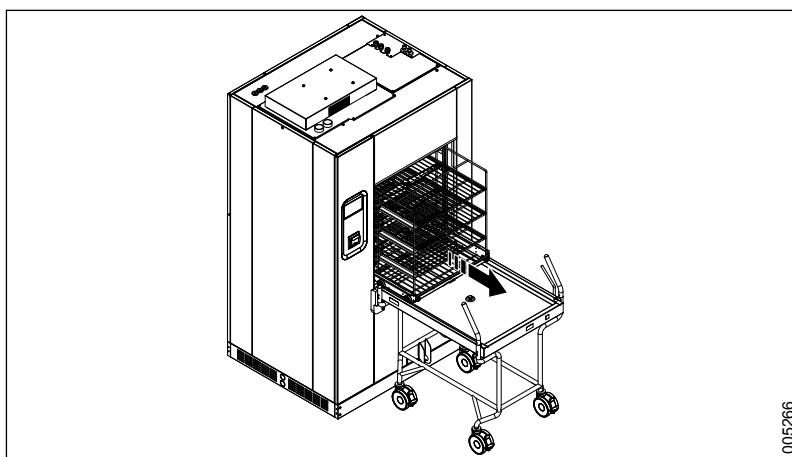

4. Undock the transport trolley.
5. Tap **Close Door** to close the door.

4.6.1 Unloading the chamber with an automatic loading system (option)

The wash cart unloads automatically onto the unloading station.

1. Dock the transport trolley to the unloading station. For more information, see section 3.2.14.2 *Transport trolley*.
2. The unloading station automatically unloads the wash cart on the transport trolley when docked. The wash cart needs to be pulled the last 10 to 15 cm manually onto the transport trolley.
3. Undock the transport trolley.

4.7 Check cleaning and disinfection results

1. Visually assess the cleanliness of the processed items, according to local requirements.
2. When processing lumen objects, verify that the instruments have remained connected during the whole process.
3. Verify, that the disinfection temperature and time were achieved by reviewing the process records according to local requirements.

4.7.1 Handling printouts (option)

If the machine is equipped with a printer, a process report is printed automatically when the program is completed. For more information, see section 6.9 *Maintenance of printer (option)*.

Remove the printout and file it in accordance with local regulations for documentation.


4.8 Holding or aborting program



If a program is interrupted and the items are not completely cleaned, the user may be exposed to contaminated items.

During a program, it may be necessary to pause the program to, for example, adjust the load.

Path: **Status View > Abort**

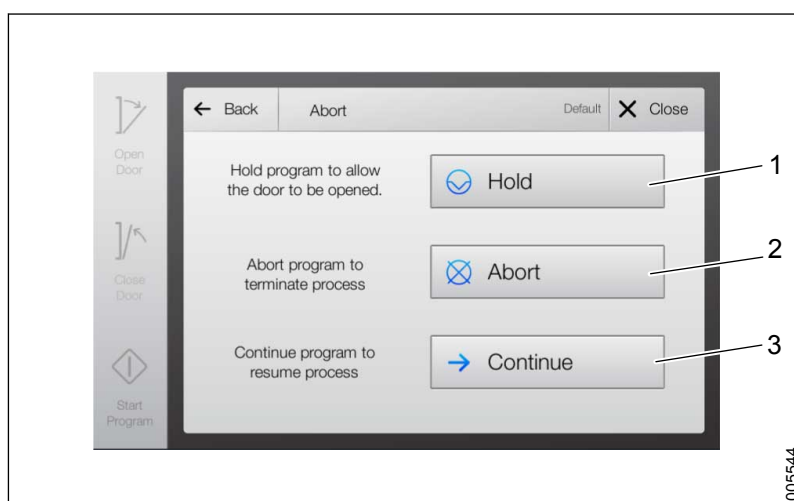

- | | |
|-----------------|--------------------|
| 1. Hold | 3. Continue |
| 2. Abort | |

Name	Description
Hold	Button for pausing the program.
Abort	Button for terminating the program.
Continue	Button for resuming the program.

Hold program (only possible in the first pre-rinse phase)

1. In the **Abort** view, tap **Hold**. The program stops and the water remains in the chamber. The door on the soiled side can be opened.
2. Tap **Continue** to resume the program from where it was stopped. The phase process time will start over again.

Open door before resuming the program

1. Tap **Open Door**.
2. Adjust the load.
3. Tap **Close Door**.

Abort program (possible at any time during a program)

1. In the **Abort** view, tap **Abort**.
Program Ends view shows until the close-down is completed. The program stops and the water is drained from the chamber. The door on the soiled side can be opened.

4.9 Change program

If the wrong program is selected, it is possible to cancel the selected program while the door is closing.

1. Tap **Open Door** while the door is closing.
2. When the door is completely open, select a new program.
3. Tap **Start program**.

4.10 Preparations for downtime

If a machine will not be used for more than 72 hours, it should be drained of all process water. Contact a service technician.

5 DETERGENTS



When handling detergents, follow the safety instructions as indicated on the detergent container label and/or in its material safety data sheet.



When selecting detergents, read and follow the instructions and recommendations of the manufacturer of the goods to be cleaned or disinfected.

The following table is a guide.

Always make decisions regarding what detergent to use for cleaning or disinfection of a specific load (material, soil and type of items) together with the detergent supplier.

Make sure the intended use for the detergent corresponds with the program and load.

Always contact a service technician for information regarding what detergent to use in a specific program phase.

Area of use	Detergent type
Thermal disinfection of non-invasive, non-critical medical devices, non-blooded waste, aids for the disabled, footwear and wash bowls.	Alkaline Enzymatic
Thermal disinfection of non-invasive, non-critical medical devices, blooded waste, container systems used to transport medical devices, aids for the disabled, footwear and wash bowls.	Alkaline Enzymatic
Thermal disinfection of invasive and critical medical devices, such as invasive equipment, soiled, rigid and non-cannulated surgical instruments.	Alkaline Enzymatic
Chemical disinfection of non-invasive, non-critical medical devices, non-blooded waste, aids for the disabled, footwear and wash bowls.	Alkaline Enzymatic
Chemical disinfection of non-invasive, non-critical medical devices, blooded waste, container systems used to transport medical devices, aids for the disabled, footwear and wash bowls.	Alkaline Enzymatic
For neutralizing alkaline residues.	Acidic

5 DETERGENTS

Area of use	Detergent type
For removal of rust, hard water scale and stains on instruments and in washer-disinfector chamber.	Acidic
For use in the final rinse to speed up the drying of the load.	Rinse aid
For use in the final rinse to lubricate moving parts of surgical instruments and prevent friction corrosion.	Lubricant

6 MAINTENANCE

This section contains checklists for maintenance intervals as well as instructions for cleaning, decontamination and re-filling detergents.

The recommended maintenance interval is mainly based on the number of programs used per day. In this manual, the specified maintenance intervals are based on an average usage rate of 10 programs per day. With a more frequent usage, the maintenance must be done more often. Other factors, such as type of load, water quality, pre-treatment and environment can also affect the need for maintenance.

Never perform service or maintenance other than described in this manual without documented experience.

6.1 Daily maintenance

Area	Check	Instructions
Exterior	Check the outer parts of the machine.	<i>6.8 Maintenance of control panel and stainless steel surfaces</i>
Detergent cabinet	Check the remaining detergent in the containers in the detergent cabinet.	<i>6.7 Maintenance of the detergent system</i>
Chamber	Check that no objects have fallen from the wash cart.	Remove loose items.
Spray arms	Check that the spray arms in the chamber rotate easily.	<i>6.4.1 Maintenance of spray arms</i>
Spray arms	Check that the holes in the spray arms are not blocked.	<i>6.4.1 Maintenance of spray arms</i>
Docking function	Check that the wash cart docks properly with the docking port.	Push in a wash cart and ensure that the docking elements align with the inlets in the wash cart manifold.
Docking function	Check that the docking springs move smoothly.	Press the docking spring in and release it. The docking spring must return to its original position.
Accessories	Check visually that the accessories have the correct function according to the manufacturer's instructions.	Visual inspection.
Accessories	Check that there are no blocked ducts, worn or missing parts that may impair the performance of the machine.	Visual inspection.
Strainer cup	Check that the strainer cup is clean and that water can flow through it.	<i>6.4.3 Maintenance of strainer</i>
Wash cart	If fine filters are used in ramps/wash carts, these should be checked and cleaned as needed.	Follow the Getinge Loading Equipment User manual.
Printer (option)	Check printout quality.	<i>3.3.7.2 Print test page (login required)</i>
Printer (option)	Check that there is sufficient paper in the printer.	<i>6.9 Maintenance of printer (option)</i>

6.2 Weekly maintenance

Area	Check	Instruction
Chamber	Check that the spray arm manifold is clean and not clogged.	<i>6.4.2 Maintenance of spray arm manifold</i>
Strainer	Check that the strainer is clean and that water can flow through it.	<i>6.4.3 Maintenance of strainer</i>
Detergent system	Check that the suction lances are intact.	Visual inspection.
	Check that the hoses are intact, flexible and soft.	Visual inspection.

6.3 Monthly maintenance


Area	Check	Instruction
Chamber	Check that the inside surfaces are free of lime scale deposits.	<i>6.5 Removing lime scale deposits in the chamber</i>

6.4 Maintenance of parts inside the chamber

6.4.1 Maintenance of spray arms




Handle the black plastic ends with care. They contain electronics for spray arm monitoring that are sensitive to mechanical influence. Blue plastic ends do not contain electronics.

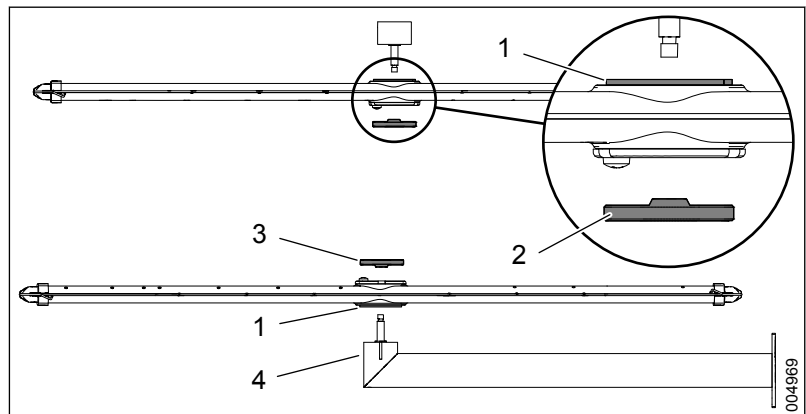

- | | |
|-------------------------|----------------------------|
| 1. Bearing ring | 3. Center nut, without lug |
| 2. Center nut, with lug | 4. Spray manifold |

The spray arm can be removed and cleaned when needed.

1. Loosen the center nut on the spray arm and remove it.
2. Remove the plastic end pieces by pressing the metal clip over the edges of the spray arm and then carefully separate the plastic end caps from the metal arm.

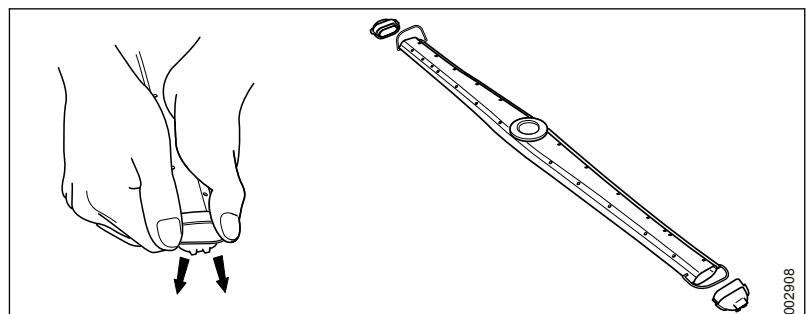

3. Using a bottle brush, clean the inside of the spray arm.

Reassemble the spray arm:

1. Reinstall the plastic end pieces on the spray arm and press the metal clips back into the groove on each side.

2. Reinstall the spray arm on the spray manifold.
The upper spray arm must be aligned with the lug of the center nut facing upwards.
Ensure that the bearing ring covers the spray manifold.
3. Check that the spray arm can rotate freely.

6.4.2 Maintenance of spray arm manifold

1. Remove the spray arm manifold from the chamber wall by loosening the nuts.
2. Remove the spray arm from the spray arm manifold; see section 6.4.1 *Maintenance of spray arms*.
3. Using a bottle brush, clean the inside of the spray arm manifold.
4. Refit the spray arm manifold.
5. Refit and check the spray arm according to the description in section 6.4.1 *Maintenance of spray arms*.

6.4.3 Maintenance of strainer

A dirty strainer could prevent water from circulating and create conditions for growth of bacteria.


WARNING!**Risk of burns**

The strainer and nearby parts may be hot.


- Beware of hot surfaces when cleaning the strainer.
 - Let the surfaces cool down before handling.
 - Wear protective gloves.
-



When removing the strainer and strainer cup, be careful not to drop loose parts under the strainer.

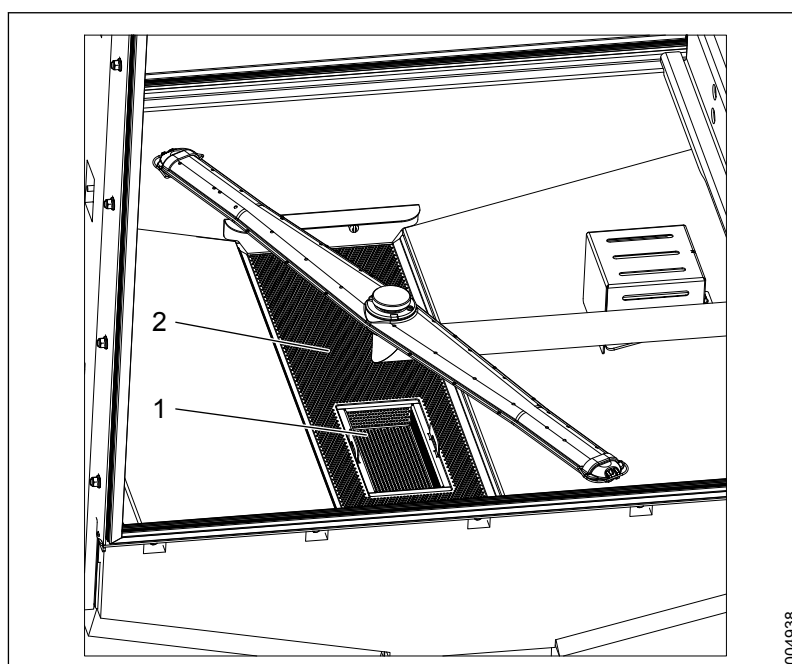

1. Strainer cup
2. Strainer

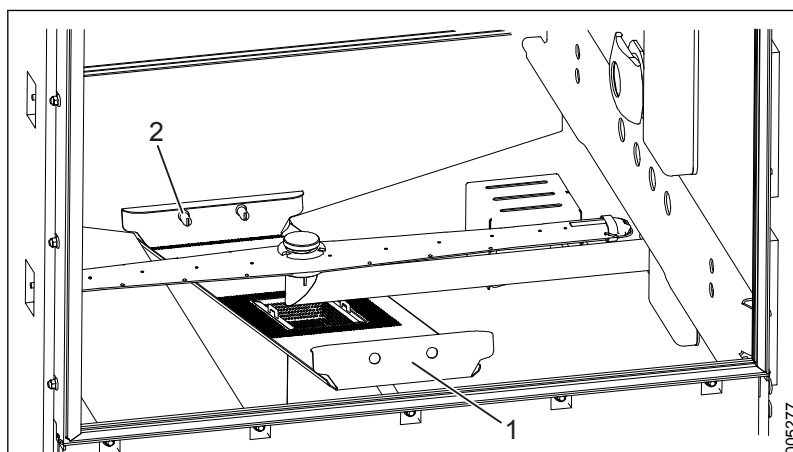

1. Strainer flange
2. Pins

Cleaning the strainer cup:

1. Check for loose items on the strainer and remove the items if necessary.
2. Brush dirt and deposits from the strainer into the strainer cup.
3. Remove and empty the strainer cup.
4. Clean the strainer cup and rinse thoroughly with clean water.
5. Put the strainer cup back in position.

Cleaning the strainer:

1. Check for loose items on the strainer and remove the items if necessary.
2. Remove the strainer by lifting it on the soiled side and pulling it outwards.
3. Clean the strainer by rinsing it with water until it is clean and water can flow through it.
4. Put the strainer back in position and make sure that the pins are fitted correctly into the holes in the strainer flange.

6.5 Removing lime scale deposits in the chamber

The program is used only to descale the machine and inserts. Do not run the **Deliming** program with items in the wash cart.

1. Pour 200 ml Getinge Clean Neutralizer Plus into a bowl and place in the chamber.
2. Load an empty wash cart into the machine.
3. Start the **Deliming** program.


6.6 Disinfecting chamber

Disinfect the chamber as follows:

1. Perform the steps in section 6.4 *Maintenance of parts inside the chamber*.
2. Run the **Self-Disinfection** program.

6.7 Maintenance of the detergent system

6.7.1 Open the detergent cabinet door

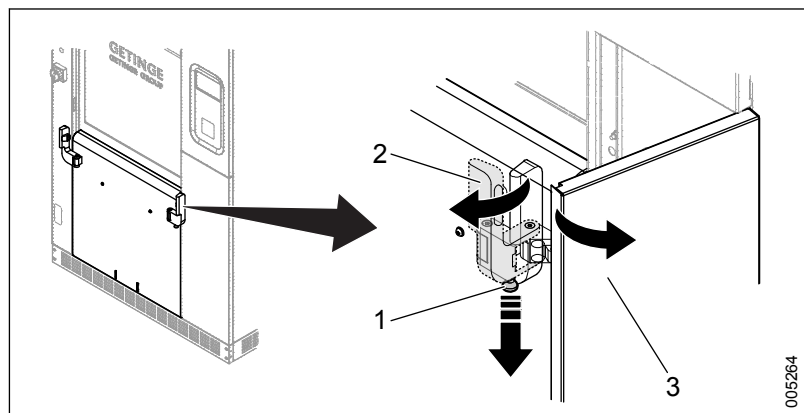

- | | |
|-----------------|---------------------------|
| 1. Locking pin | 3. Detergent cabinet door |
| 2. Docking part | |

1. Pull the locking pin downwards.
2. Rotate the docking part outwards.
3. Open the cabinet door.

6.7.2 Check for damage

Visually inspect the hoses and the suction lances:

- Ensure that there is no leakage.
- Ensure that the hoses are not cracked or damaged.
- Ensure that the hoses are flexible and soft.

If any of the above is not fulfilled, contact service personnel.

6.7.3 Refill detergent



**Keep the suction lances and hoses clean.
Remove excess and/or crystallized detergent.**



When changing from one type of detergent to another, the dosing pump must be recalibrated.



If a detergent label or dosage tag is missing or if the product name on the label is not the same as the product name on the detergent container, contact service personnel.

1. Remove the empty container from the detergent cabinet. For more information, see section 6.7.1 *Open the detergent cabinet door*.
2. Remove the suction lance from the empty container.
3. Clean the suction lance with a cloth. Ensure that it is free of crystallized liquid.
4. Insert the suction lance into a full container. Ensure that the container contains the correct detergent. For more information, see section 3.2.3.1 *How detergents are labeled*.
5. Place the full container in the detergent cabinet.
6. Close the door.

6.8 Maintenance of control panel and stainless steel surfaces

Control panel

1. Wipe off stains and dirt with a wet or dry cloth and an alcohol-based cleaning agent or rubbing alcohol.

Stainless steel surfaces


1. Wipe off stains and dirt using a wet or dry cloth and detergent for stainless steel.

6.9 Maintenance of printer (option)

The built-in thermal printer (option) uses heat-sensitive paper for printing. No ink cartridge needs to be replaced.



The control system alerts with a warning on the touchscreen when the printer is low on paper.

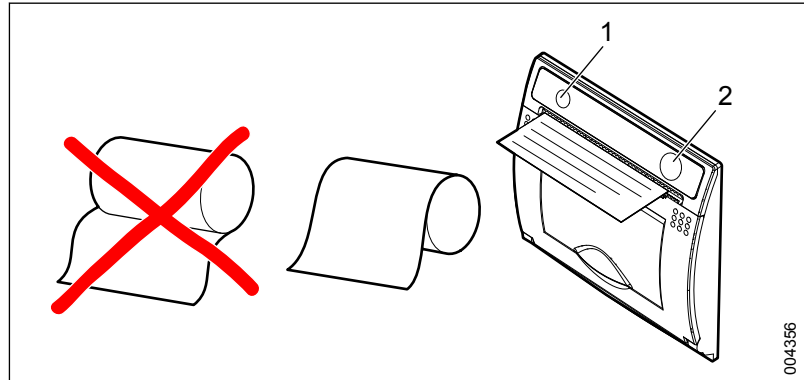

- 1. Indication lamp
- 2. Feed paper button

1. Open the lid on the printer.
2. If the back side of the paper has red markings, the paper is about to run out and needs to be replaced.
3. Remove the existing roll and replace it with a new one. Position the roll in the correct direction.
4. Close the lid.
5. Feed paper forward until blank paper appears.

7 TROUBLESHOOTING

This section contains information and instructions for troubleshooting a malfunctioning unit.

7.1 System alerts

7.1.1 Handle alarms (login required)

Alarm during ongoing program

1. Tap **Mute Alarm** to shut down the alarm signal.
2. Note the displayed message, error code and phase. Tap **Guide** for more details.
3. Take action according to the Alarms and warnings list. For more information, see section *7.1.3 Alarms and warnings list*.
4. Login and tap **Confirm Alarm**.
5. Ensure that the doors are closed. If they are not closed, tap **Close Door**.
6. Tap **End Program** and **Program Ends** is displayed while the machine is drained of water. When the program shut-down is completed, the door on the soiled side can be opened and **Goods not processed** is displayed on the touch screen.
7. Because the load is not processed, a new program must be selected. Select a new program to process the load correctly, tap **Open Door** in the **Goods not processed** view and the **Select Program** view appears.
8. Select an appropriate program and tap **Start program**.

Alarm during standby

1. Follow steps 1-4 from Alarm during ongoing program.
2. Tap continue and **Select Program** view appears.

7.1.2 Handle warnings

A warning indicates a minor problem. The program is not stopped.

1. Note the displayed message, error code and phase. Tap **Guide** for more details.
2. Tap **Confirm**.
3. Take action according to the Alarms and warnings list; see section *7.1.3 Alarms and warnings list*.

7.1.3 Alarms and warnings list



If the action in the following list does not resolve the problem, contact service personnel.

* red = alarm, yellow = warning				
Code	*	Displayed message	Cause	Action
11002		Power failure	The program has been aborted due to power failure.	Confirm the alarm and restart the program.
11003		I/O module fault	I/O module error.	Contact service personnel.
11004		Process stop	The process stop has been activated.	Reset the process stop and confirm the alarm.
11010		SS door lock	The door on the soiled side does not lock.	Ensure that nothing is blocking the door.
11011		SS door unlock	The door on the soiled side does not unlock.	Ensure that nothing is blocking the door.
11012		SS door close	The door on the soiled side does not close.	Ensure that nothing is blocking the door.
11013		SS door open	The door on the soiled side does not open.	Ensure that nothing is blocking the door.
11015		CS door lock	The door on the clean side does not lock.	Ensure that nothing is blocking the door.
11016		CS door unlock	The door on the clean side does not unlock.	Ensure that nothing is blocking the door.
11017		CS door close	The door on the clean side does not close.	Ensure that nothing is blocking the door.
11018		CS door open	The door on the clean side does not open.	Ensure that nothing is blocking the door.
11020		Slow fill chamber	The filling of the chamber takes too long time.	Contact service personnel.
11021		Slow drain chamber	The draining of the chamber takes too long time.	Contact service personnel.
11022		Low level chamber	The water level in the chamber is too low.	Contact service personnel.
11023		Low circ. press	The pressure in the circulation pump is too low.	Ensure that the wash cart is correctly docked in the chamber.
11024		High circ. press	The pressure in the circulation pump is too high.	Ensure that the spray arms, nozzles or spray pipes are not blocked by objects or clogged. For more information, see section 6 MAINTENANCE.

* red = alarm, yellow = warning				
Code	*	Displayed message	Cause	Action
11025		Water leakage dryer	Water has been detected in the drying system.	Contact service personnel.
11026		Slow drain D. tank	The draining of the drain tank is taking too long. Only applicable for S-8668T.	Contact service personnel.
11027		Slow spray wing	One or more of the spray arms are rotating too slowly. If the problem is not solved within 15 minutes after the warning, an alarm is triggered.	Ensure that the spray arms, nozzles or spray pipes are not blocked by objects or clogged. For more information, see section 6 <i>MAINTENANCE</i> .
11028		Circ. pump fail	Error in the circulation pump.	Contact service personnel.
11029		Drain valve fail	The drain valve is open. Only applicable for S-8668T.	Contact service personnel.
11030		Low temp chamber	The chamber temperature is too low.	Contact service personnel.
11031		High temp chamber	The chamber temperature is too high.	Contact service personnel.
11032		Temp deviation	Faulty temperature sensor.	Contact service personnel.
11033		Disinfection fail	The disinfection phase was not successful.	Contact service personnel.
11034		Slow fill proc. tank	The filling of the process tank is taking too long. Only applicable for S-8668T.	Contact service personnel.
11035		Slow drain proc. tank	The draining of the process tank is taking too long. Only applicable for S-8668T.	Contact service personnel.
11036		Proc. tank sensor fail	Faulty process tank level sensor. Only applicable for S-8668T.	Contact service personnel.
11037		High conductivity	Process water has too high conductivity reading.	Contact service personnel.
11040		Dosage pump 1	Dosing failure pump 1.	Ensure that the level of detergent is sufficient.
11041		Dosage pump 2	Dosing failure pump 2.	Ensure that the level of detergent is sufficient.
11042		Dosage pump 3	Dosing failure pump 3.	Ensure that the level of detergent is sufficient.
11043		Dosage pump 4	Dosing failure pump 4.	Ensure that the level of detergent is sufficient.
11050		Slow fill booster	The filling of the booster tank is taking too long.	Contact service personnel.

* red = alarm, yellow = warning				
Code	*	Displayed message	Cause	Action
11051		Slow drain booster	The draining of the booster tank is taking too long.	Contact service personnel.
11052		Low temp booster	The booster tank temperature is too low.	Contact service personnel.
11053		High temp booster	The booster tank temperature is too high.	Contact service personnel.
11060		Low temp dryer	The drying unit temperature is too low.	Contact service personnel.
11061		High temp dryer	The drying unit temperature is too high.	Contact service personnel.
11062		Low press dryer	The drying unit pressure is too low.	Contact service personnel.
11063		High press dryer	The drying unit pressure is too high.	Contact service personnel.
11070		S Disinfection fail	Disinfection fail alarm.	Contact service personnel.
11071		S Temp deviation	Temperature deviation alarm.	Contact service personnel.
11072		S Circ. press	Circulation pump pressure alarm.	Contact service personnel.
10073		S H conductivity	High conductivity alarm.	Contact service personnel.
11074		S Dosage	Dosing fail alarm.	Contact service personnel.
11075		Supervisor error	Supervisor error alarm.	Contact service personnel.
11080		AI fail-Chamb temp	Sensor failure.	Contact service personnel.
11081		AI fail-Dryer temp	Sensor failure.	Contact service personnel.
11082		AI fail-Chamb I temp	Sensor failure.	Contact service personnel.
11083		AI fail-Boost temp	Sensor failure.	Contact service personnel.
11084		AI fail-Circ press	Sensor failure.	Contact service personnel.
11085		AI fail-Humidity	Sensor failure.	Contact service personnel.
11086		AI fail-Conductivity	Sensor failure.	Contact service personnel.
11087		AI fail-Dryer press	Sensor failure.	Contact service personnel.
13001		Chemical 1 low level	The level in the detergent container is too low.	Refill detergent or replace the container.
13002		Chemical 2 low level	The level in the detergent container is too low.	Refill detergent or replace the container.
13003		Chemical 3 low level	The level in the detergent container is too low.	Refill detergent or replace the container.

* red = alarm, yellow = warning				
Code	*	Displayed message	Cause	Action
13004		Chemical 4 low level	The level in the detergent container is too low.	Refill detergent or replace the container.
13007		Scanner read error	The machine's scanner was not able to read the bar code on the wash cart.	Pull out the wash cart and ensure that the barcode is clean and not worn. Push the wash cart back in and start the program. If the problem persists, use another wash cart.
13008		ID tag not detected	ID tag on wash cart cannot be found. The tag is used to indicate how many spray arms must be monitored for rotation.	Pull the cart out of the chamber. Make sure the ID tag is in place. Push the cart in again and start the program.
13010		Slow spray wing	One or more of the spray arms are rotating too slowly.	Open the door and check the wash cart. Ensure that the spray arms, nozzles or spray pipes are not blocked by objects. Close the door and continue the program. (The phase time will restart).
13021		Chemical 1 empty	The detergent container is empty.	Refill detergent or replace the container.
13022		Chemical 2 empty	The detergent container is empty.	Refill detergent or replace the container.
13023		Chemical 3 empty	The detergent container is empty.	Refill detergent or replace the container.
13024		Chemical 4 empty	The detergent container is empty.	Refill detergent or replace the container.
93000		Config. File: Read Error	Control system error.	Contact service personnel.
93001		System Time Monitor Err	Control system error.	Restart the machine.
93002		Data Integrity Error	Control system error.	Restart the machine.
93003		Internal System Error	Control system error.	Restart the machine.
93004		Program Database Error	Control system error.	Restart the machine.
93005		File System: Mem Low	Control system memory low.	Contact service personnel.
93006		File System: Mem Full	Control system error.	Contact service personnel.
93007		Program Import Fail	Import rejected.	Contact service personnel.
93008		Supervisor Version Error	Control system error.	Contact service personnel.
93100		Users File: Read Error	Control system error.	Contact service personnel.


* red = alarm, yellow = warning				
Code	*	Displayed message	Cause	Action
93101		Rights File: Read Error	Control system error.	Contact service personnel.
93102		X20Config File: Read Err	Control system error.	Contact service personnel.
93103		License File: Error	Failed to read license file for software addition during the boot process.	Contact service personnel.
93104		License File: Error	The panel serial number and the serial number in the license file do not match. The license file has been copied/reset from another panel.	Contact service personnel.
93200		Process Report Error	System not able to print a process report.	Print the report again. For more information, see section 3.3.7.1 <i>Print emergency print-outs (login required)</i> .
93201		Netw. Printer: Conn Err	Connection problems with the network printer.	Contact local IT-support to check the network connection.
93202		Fascia Printer: Conn Err	Connection problems with the fascia printer.	Contact service personnel.
93203		USB Storage: Conn Err	Connection problems with the USB device.	Contact service personnel.
93204		Netw. Storage: Conn Err	Connection problems with the network storage device.	Contact local IT-support to check the network connection.
93206		USB Storage: Write Err	The report is not saved on the USB stick	Change the USB memory stick.
93207		USB Storage: Mem Full	The USB stick does not have sufficient space for storage.	Change the USB memory stick.
93208		USB Storage: Mem Low	The USB stick does not have sufficient space for storage.	Change the USB memory stick.
93209		Netw. Storage: Write Err	The report has been archived and a new attempt to send the report will be made as soon as a connection is re-established.	Contact local IT-support to check the network connection.
93210		Fascia Printer: Write Err	The control system is unable to print the process report on the fascia printer.	Check the printer paper roll. Refill paper if needed. Ensure that the printer lid is closed. Print the report again.

* red = alarm, yellow = warning				
Code	*	Displayed message	Cause	Action
93211		Fascia Printer: No Paper	No paper in the Fascia printer.	Check the printer paper roll. Refill paper if needed. See section 6.9 <i>Maintenance of printer (option)</i> . Ensure that the printer lid is closed. Print the report again. For more information, see section 3.3.7.1 <i>Print emergency print-outs (login required)</i> .
93212		Netw. Printer: Write Err	The report has been archived and a new attempt to send the report will be made as soon as a connection is re-established.	Check the printer paper roll. Refill paper if needed. Check the printer toner level. Contact local IT-support to check the network connection.
93214		Netw. Printer: Report OK	Process report is successfully sent to the network printer.	No action required..
93215		Fascia Printer: Report OK	Process report is successfully sent to the fascia printer.	No action required..
93216		USB Storage: Report OK	Process report is successfully sent to the USB device.	No action required..
93217		Netw. Storage: Report OK	Process report is successfully sent to the network storage device.	No action required..
93301		Barcode Reader: Conn Err	Failed to connect to the bar code reader via USB.	Contact service personnel.
93508		T-DOC: System Unavailable	T-DOC error	Contact local IT-support to check the network connection.
93509		T-DOC: Connection Lost	T-DOC error	Contact local IT-support to check the network connection.
93551		GO: System Unavailable	Getinge Online connection error	Contact local IT-support to check the network connection.
93552		GO: Connection Lost	Getinge Online connection error	Contact local IT-support to check the network connection.
93553		GO: Hostname Fail	Getinge Online connection error	Contact local IT-support to check the network connection.

* red = alarm, yellow = warning				
Code	*	Displayed message	Cause	Action
93555		GO:Proxy Login Fail	Getinge Online connection error	Contact local IT-support to check the network connection.
93556		GO:Gateway Error	Getinge Online connection error	Contact local IT-support to check the network connection.
93557		GO:Gateway Config	Getinge Online connection error	Contact local IT-support to check the network connection.

7.1.4 Handle information message

An information message displays a condition that must be resolved before a program can be started.

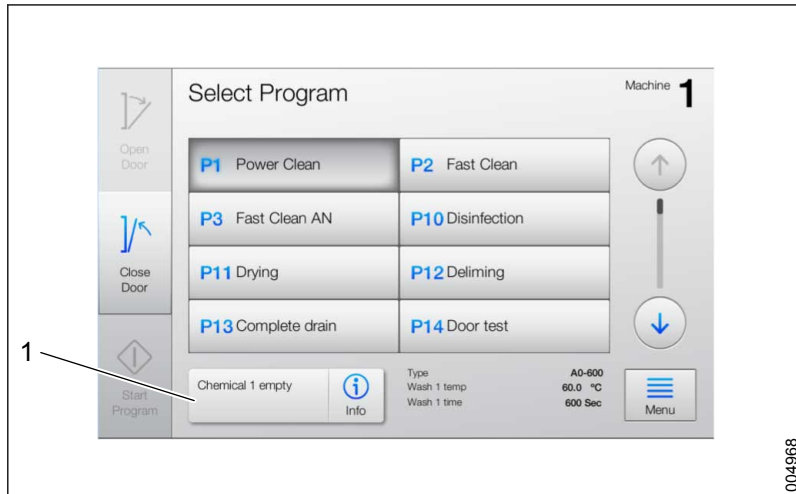

1. Information message

1. Check the information message.
2. Tap **Info** to get more detailed information.
3. Resolve the condition.
4. Start the program. No confirmation is required.



If the action in the following list does not resolve the problem, contact service personnel.

Displayed message/text	Cause	Action
CS Door not locked	The door on the clean side is not closed and locked.	Close the clean side door.
Scanner read error	The machine's scanner was not able to read the barcode on the wash cart.	Rescan the barcode or select a program on the touchscreen.
ID tag not detected	ID tag on wash cart cannot be found. The tag is used to indicate how many spray arms must be monitored for rotation.	Pull out the wash cart and check the ID tag. Push the wash cart back in and start the program. If the problem persists, use another wash cart.
Alarm active	An alarm is triggered.	For more information, go to Menu > Alarms & Messages .
No program selected.	No program selection made.	Select a program on the touchscreen.
Config error	Disinfection is not correctly configured.	Contact service personnel.
Config error	Incorrect water selection configured. Only applicable for S-8668T.	Contact service personnel.
Log not ready	A log-file is being prepared.	Wait until the task is completed.
Config error	The dosing is not correctly configured.	Contact service personnel.
Chemical 1 empty	The detergent container is empty.	Refill detergent or replace the container.
Chemical 2 empty	The detergent container is empty.	Refill detergent or replace the container.
Chemical 3 empty	The detergent container is empty.	Refill detergent or replace the container.
Chemical 4 empty	The detergent container is empty.	Refill detergent or replace the container.
Manual output active	One or more outputs are active.	Contact service personnel.

7.2 General troubleshooting



If the action in the following list does not resolve the problem, contact service personnel.

Incident	Cause	Action
The wash cart does not dock properly.	The wash cart is not placed correctly in the chamber.	Place the wash cart according to instructions. For more information, see section 4.4 <i>Loading the chamber</i> .
The cleaning or drying result is not satisfactory.	The load is not placed correctly on the wash cart.	Reload the items according to instructions. For more information, see section 4.4 <i>Loading the chamber</i> . Clean again.
The cleaning result is not satisfactory.	The load is too soiled for the selected program.	Prepare the load according to instructions. For more information, see section 4.3 <i>Preparations</i> . Select an appropriate program according to the level of contamination.
The cleaning result is not satisfactory.	The selected detergent is not correct.	Contact service personnel.
The cleaning result is not satisfactory.	The selected accessories are not working correctly.	Use only approved accessories and follow manufacturers instructions. For more information, see section 3.4 <i>Accessories</i> .
The drying result is not satisfactory.	No rinse agent is used.	Contact service personnel. Add rinse agent to the program.
The drying result is not satisfactory.	The drying time is too short.	Contact service personnel. Extend the drying time. This requires a login.
The cleaning result is not satisfactory.	Cleaning process parameters need optimization.	Contact service personnel.
The cleaning result is not satisfactory.	The items/instruments are too worn out to be cleaned effectively.	Replace the items/instruments.
The printer light is flashing.	Paper error.	The paper must be corrected or the printer refilled. For more information, see section 6.9 <i>Maintenance of printer (option)</i> .
No printout comes out and no warning is triggered.	Printer is out of paper.	Refill printer. For more information, see section 6.9 <i>Maintenance of printer (option)</i> and 3.3.7.1 <i>Print emergency print-outs (login required)</i> .
The washer-disinfector does not start.	The main power is not switched on.	Switch on the main power on the machine.
The washer-disinfector does not start.	The isolator switch is not switched on.	Switch on the isolator switch.
The washer-disinfector does not start.	A fuse has been tripped.	Contact service personnel.

Incident	Cause	Action
The washer-disinfector does not start.	The residual current device has been tripped.	Contact service personnel.
The washer-disinfector does not empty the water at the end of a program.	F10 UH02 has been tripped.	Contact service personnel.
The washer-disinfector acts in an unexpected way (noise, smell etc.).	An internal component has failed.	Contact service personnel.
The washer-disinfector acts in an unexpected way (noise, smell etc.).	The connections to the washer-disinfector might not be correct according to specifications.	Close off all connections to the machine (water, electricity, steam etcetera). Contact service personnel.
The touchscreen is black.	The screen saver is enabled; the touchscreen is in idle mode.	Tap on the touchscreen to light it up.
The touchscreen is black.	The main power switch is not switched on.	Switch on the main power on the machine.
The touchscreen is black.	The isolator switch is not switched on.	Switch on the isolator switch.
The touchscreen is black.	A fuse has been tripped.	Contact a local electrician.
The touchscreen is black.	The touchscreen is in start-up mode. During the start-up procedure, the touchscreen is black for approx. 15 seconds.	Wait for the touchscreen to start up.
The touchscreen displays a warning that the USB memory is almost full.	The USB memory is almost full.	Transfer the data to another storage device.
The scanner was not able to read the barcode on the wash cart.	The barcode is dirty or worn out.	Clean the barcode and rescan it. If the barcode can't be read, use the touchscreen to select a program.
Need to change the selected program.	An incorrect program is selected.	For more information, see section 4.9 <i>Change program</i> .
Need to hold or abort a program.	The load needs to be adjusted.	For more information, see section 4.8 <i>Holding or aborting program</i> .

7.3 Washing efficiency

If the washing efficiency deteriorates over time, follow local regulations or contact service personnel.


7.4 Drying efficiency

If the drying efficiency deteriorates over time, follow local regulations or contact service personnel.

7.5 Restart the machine

To restart the machine, use the main power switch to do the following:

1. Turn the main power off (0).
2. Wait until the touchscreen is deactivated.
3. Turn the main power on (1).

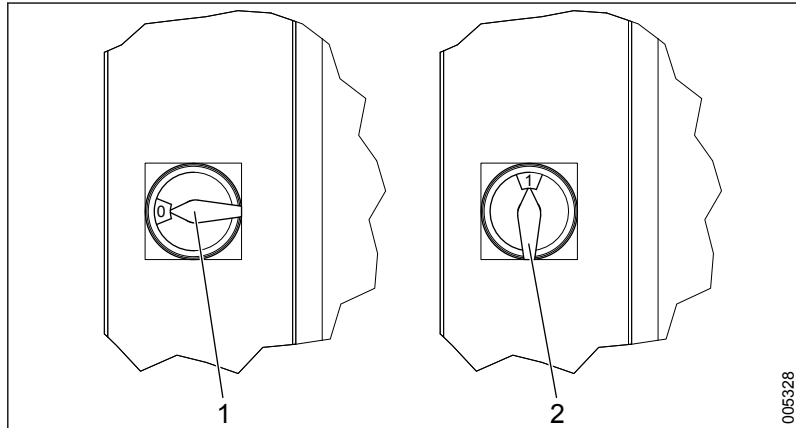

1. Main power off
2. Main power on

8 HANDLING OF WORN-OUT PRODUCTS

This section describes how to dispose of worn out products and parts.

8.1 WEEE directive

- This product and its accessories comply with the requirements in the WEEE directive.
- WEEE stands for "Waste Electrical & Electronic Equipment".

8.2 Disposal of detergent or detergent containers

- Left-over detergent and used detergent containers must be disposed of at recycling stations in accordance with national laws and local regulations.
- If required, contact your detergent supplier for more information.

8.3 Disposal of the machine, accessories or spare parts



Before recycling, make sure that the machine is disinfected.

- Worn-out products must be disposed of at recycling stations in accordance with national laws and local regulations.
- If required, contact your reseller for more information.

8.3.1 List of materials

Each type of material must be handled and recycled according to current national laws on waste handling. There is no mercury in the machine or its components.

Materials	Approximate share of materials					
	S-8666		S-8668		S-8668T	
	kg	%	kg	%	kg	%
Stainless steel	299.0	70	313.0	70	342.3	70
Steel	7.9	2	8.2	2	9.0	2
Cast iron	2.8	< 1	2,9	1	3.2	< 1
Copper	15.3	4	16.0	4	17.5	4
Plastic	69.8	16	73.1	16	80.0	16
Glass	18.4	4	19.3	4	21.1	4
Rubber/Silicone	0.3	< 1	0.3	< 1	0.4	< 1
Wiring *	5.3	1	5.6	1	6.1	1
Circuit boards *	2.5	< 1	2.7	< 1	2.9	< 1
Other components **	7.0	2	7.3	1	8.0	2
Total, all materials	428.3	100	448.4	100	490.5	100
* contains environmentally hazardous substances such as lead, cadmium, PVC and flame retardants and must be recycled as electronic scrap. ** consists of pumps, door motors and other parts that cannot be weighed separately.						

Packing material

The machine is delivered on wooden pallets and the machine parts are covered with corrugated cardboard or wrapped in plastic.

9 ADDRESSES

Australia

Getinge Australia PTY Ltd
PO Box 50
Bulimba QLD 4171
Australia
Cust Support: 1300 155 500
Phone: +61 7 3399 3311
Fax: +61 7 3395 6712
E-mail:
info@getinge.com.au

China

Getinge (Suzhou) Co.Ltd
No.158, Fang Zhou Road,
Suzhou Industrial Park
215 021 Suzhou
Jiangsu Province P.R. China
Phone: +86-51 262 838 966
E-mail:
info@getinge.com.cn

Finland

Getinge Finland AB
Niittykatu 8
FI-02200 Espoo
Finland
Phone: +35-89 6824 120
E-mail:
getinge.finland@getinge.com

Italy

Getinge S.p.A
via Poggio Verde, 34
00148 Rome
Italy
Phone: +39-06 656 631
E-mail:
info@getinge-it

Belgium

Getinge NV
Vosveld 4 B-2
B-2110 Wijnegem
Belgium
Phone: +32-33 542 865
E-mail:
info@getinge.be

China

Getinge Shanghai Trading
Co.Ltd.
No.3, Lane 128, Lin Hong Road
Changing District
200335 Shanghai China
Phone: +86-21 6197 3999
E-mail:
info@getinge.com.cn

France

Getinge France SAS
BP 49, avenue du Canada
ZA de Courtaboeuf
Les Ulis, FR-91942
France
Phone: +33-1 64 86 89 00
E-mail:
getinge.france@getinge.fr

Japan

Getinge Group Japan K.K.
SPHERE TOWER TENNOZ
23rd Fl.
2-2-8 Higashi Shinagawa
Shinagawa, Tokyo 1400002
JAPAN
Phone: +81-3- 6863-8004

Canada

Getinge Canada Ltd
6685, Millcreek Drive
Unit 3-5
Mississauga, Ontario L5N5M5
Canada
Phone: +1-905-629-8777
E-mail:
customercare@getinge.ca

Denmark

Getinge Danmark A/S
Industriparken 44 B
DK-2750 Ballerup
Denmark
Phone: +45 45 93 27 27
Fax: +45 45 93 41 20
E-mail:
getinge.denmark@getinge.com

Germany

Getinge Vertrieb & Service
GmbH
Kehler Strasse 31
764 37 RASTATT
Germany
Phone: +49-7222 932 306
Fax: +49-7222 932 597
E-mail:
info.inco-de@getinge.com

Netherlands

Getinge b.v.
Biezenwei 21
4004 MB Tiel
The Netherlands
Phone: +31 (0)344-809900
Fax: +31 (0)344-640885
E-mail:
info@getinge.nl

Norway

Getinge Norge A/S
Strandveien 13
1366 Lysaker
Norway
Phone: +47-23 03 52 00
E-mail:
info@getinge.no

Spain

Getinge Iberica SL
C/ marie Curie, 5
Edificio Alfa, Planta 6a
28521 Rivas Vaciamadrid
Spain
Phone: +34-916 78 26 26
E-mail:
administracion@getinge.es

United Kingdom

Getinge UK Ltd
1 Pembroke Avenue,
Waterbeach
Cambridge, CB25 9QP
United Kingdom
Phone: +44 (0) 1223 861665
E-mail:
sales@getinge.co.uk

Poland

Getinge Poland sp. z o. o.
ul. Osmańska 14,
PL-02-823 Warsaw
Poland
Phone: +48-22 882 06 26
E-mail:
office@getinge.pl

Sweden

Getinge Sverige AB
P O Box 69
SE-310 44 Getinge
Sweden
Phone: +46 10 335 30 00
E-mail:
info@getinge.com

USA

Getinge USA Inc.
1777 East Henrietta Road
Rochester, NY 14623-3133
USA
Phone: +1-800-475-9040
E-mail:
info.northamerica@getinge.com

Singapore

Getinge International Far East
Pte. Ltd
20 Bendemeer Road,
#06-02, Cyberhub Building
Singapore, SG-339914
Phone: +65-6396 7298
E-mail:
info.sg@getinge.com

Switzerland

Getinge ALFA AG
Weidenweg 17
4310 Rheinfelden
Switzerland
www.getingealfa.ch
Phone: +41-61 836 15 15
E-mail:
info@alfo.ag

Alphabetical index

About this manual	7
Access rights	28
Accessories	40
Accident or incident reporting	14
ADDRESSES	75
Alarms	40
Alarms and warnings list	61
Ambient illumination	40
Automatic loading system (option)	27
Barcode scanner (option)	24
Basic views	32
Before using the machine	41
Calibrate touchscreen (login required)	36
Carts	26
Chamber	25
Change program	48
Change settings	35
Check cleaning and disinfection results	46
Check for damage	57
Configurations	15
Control panel	17
Control system	27
Daily maintenance	51
Declaration of Conformity	9
DESCRIPTION	15
Detergent cabinet	21
DETERGENTS	49
Disclaimer	9
Disinfecting chamber	57
Disposal of detergent or detergent containers	73
Disposal of the machine, accessories or spare parts	73
Docking	25
Door	24
Door pinch protection	12
Drying efficiency	71
Electrical outlet (option)	25
Electromagnetic environment	40
General troubleshooting	70
General warnings	13
Guide rail	25
Guides for external loading equipment (option)	24
Guides for transport trolley	23
Handle alarms (login required)	60
Handle information message	68

Handle warnings	60
HANDLING OF WORN-OUT PRODUCTS	73
Handling printouts (option)	46
Holding or aborting program	47
How detergents are labeled	22
Information text	38
Intended application areas	8
Intended use statement	7
Intended user	8
Intended user environment	8
INTRODUCTION	7
Isolator switch	12
Keyboard	29
List of materials	74
Loading the chamber	43
Loading the chamber with an automatic loading system (option)	43
Loading the wash cart	42
Lower guide for transport trolley	23
Main power switch	24
MAINTENANCE	51
Maintenance of control panel and stainless steel surfaces	58
Maintenance of parts inside the chamber	53
Maintenance of printer (option)	59
Maintenance of spray arm manifold	54
Maintenance of spray arms	53
Maintenance of strainer	55
Maintenance of the detergent system	57
Manufacturer	9
Menu	33
Menu tree	31
Monthly maintenance	52
Navigation and commands	29
Open the detergent cabinet door	57
OPERATION	41
Outlet for external loading equipment (option)	24
Password protection	12
Personal protection	13
Potential misuse	9
Preparations	41
Preparations for downtime	48
Print emergency print-outs (login required)	37
Print test page (login required)	37
Printer (option)	18

Printing	37
Process stop (option)	24
Product description	16
Product labels	10
Product liability	9
References	11
Refill detergent	58
Removing lime scale deposits in the chamber	56
Restart the machine	72
SAFETY	12
Safety functions	12
Safety labels on the machine	13
Safety statements	12
Select Program view (start view)	32
Service life	11
Set screensaver (login required)	35
Set status view	35
Set touchscreen and print-outs language (login required)	36
Set touchscreen brightness (login required)	35
Spray arm	25
Spray arm manifold	26
Starting program	44
Starting program with a barcode scanner (option)	44
Starting program with an automatic loading system (option)	45
Status view (during program)	34
Strainer	26
Symbols in this manual	7
System alerts	37, 60
Terms used in this manual	7
Touchscreen	17
Touchscreen commands	30
Transport trolley	27
TROUBLESHOOTING	60
Unloading the chamber	45
Unloading the chamber with an automatic loading system (option)	46
USB (option)	19
User Manual	1
Warnings	39
Warranty	9
Wash cart	26
Washing efficiency	71
WEEE directive	73
Weekly maintenance	52
View alarm and message information	37
View system and machine information	37
View system information	37
Workflow	41



Manufacturer:

Getinge Disinfection AB
Ljungadalsgatan 11
352 46 Växjö
SWEDEN

www.getinge.com

2021-06-07 EN REV.H